

PH-ABT-HC-RFC7S

Images



Product Description

The Premier Pharmacy Dual Temp Combination Refrigerator/Freezers feature separate microprocessor temperature controller allowing for precise temperature control, verification and recovery. An interior fan driven, forced-draft air circulation system provides excellent temperature uniformity and recovery after door openings
REFRIGERATORS ONLY.

Critical samples and supplies are security backed with an extensive package of audible and visual alarms.

*Please note this unit requires 2 outlets.

Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure. Periodic defrosting of the freezer is required.

Certifications



General Description and Application

Storage capacity (cu. ft)	5.2 cu. ft. capacity (Refrigerator) / 1.7 cu. ft. capacity (Freezer)
Door	2 exterior swing doors (solid door refrigerator/solid door freezer), right hinged, not reversible
Shelves	3 adjustable shelves refrigerator/2 adjustable shelves freezer
Freezer Compartments	1
Baskets	Non-applicable
Mounting	Leveling Legs and stacking kit
Interior lighting	LED interior light refrigerator
Airflow management	Forced draft air circulation, refrigerator
External probe access	Two probe access ports (3/4"); one, refrigerator/one, freezer
Insulation	High density urethane foam cabinet insulation
Exterior materials	Steel
Access Control	Keyed door lock
General warranty	Two (2) years parts and labor warranty
Compressor Parts Warranty	Five (5) years compressor warranty
Product Weight (lbs)	135
Shipping Weight (lbs)	215
Rated Amperage	1.3 Amps refrigerator, 1 Amp freezer
Power Plug/Power Cord	NEMA 5-15 Plug, conforms to UL471 requirements
Facility Electrical Requirement	110-120V AC: 15 A (minimum)
Agency Listing and Certification	UL-C-UL Listed

Refrigeration System

Refrigerant	Hfc-free refrigerant (r600a)
Defrost	Manual defrost freezer / cycle defrost refrigerator

Controller, Configuration, Alarms and Monitoring

Controller technology	Microprocess temperature controller
Calibration	Calibrated using a NIST traceable device, certificate included
Adjustable Temperature Range	Refrigerator: 36°F – 46°F (2°C -8°C); Freezer: 5°F – -13°F (-15°C - -25°C)
External alarm connection	Remote alarm contacts
Alarms	Audible and visual high and low temperature alarms

Other	Pyxis®, Omnicell® and AcuDose RX® compatible
Warranty Disclaimer	Indoor use only, +18C to +26C (+65F to +78F), <70% RH.

Specifications are subject to change.


Requires 2 outlets.

Models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure.

Periodic defrosting of the freezer is required.

Dimensions

	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	23.625"	24"	53"	22"	46"
Interior	17.75" upper / 19.875" lower	15.75" upper / 16" lower	13.75" upper / 26" lower		

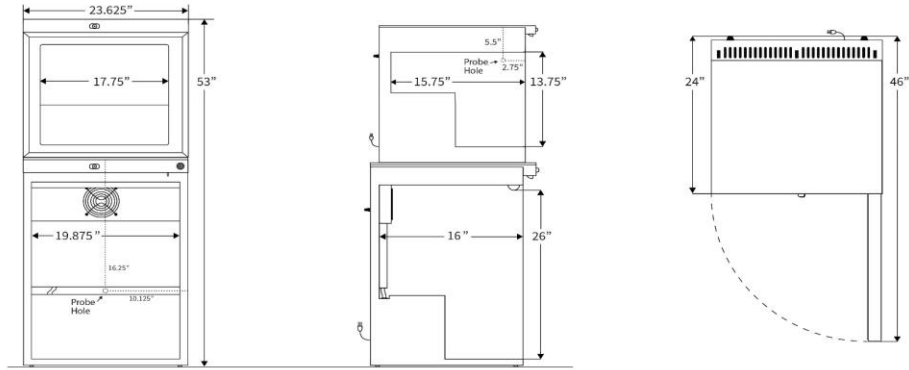

Contact

Customer Service
Technical Service

800-648-4041 Option 3
800-648-4041 Option 5, Parts Option 4

customerservice@horizontalscientific.com
technicalservice@horizontalscientific.com