

# ACR1012GLHD



53.13" x 23.38" x 28.33" (H x W x D)

Upright healthcare all-refrigerator with self-closing glass door, front lock, digital controls, and external temperature readout

## Highlights:

External temperature readout of the current and high/low temperatures

User-friendly front control panel keeps the thermostat and alarm in reach without having to open the door

Self-closing door protects interior products from ambient temperatures


## Product Features:

|                              |                                                                                                                                                                                                                                                 |
|------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| EQTemp Medical Refrigerator  | All-refrigerator designed to operate between 2 and 8°C for a variety of storage application in medical offices and other healthcare facilities                                                                                                  |
| Temperature Range: 2° to 8°C | Preset to 5°C as recommended for many storage applications                                                                                                                                                                                      |
| External temperature readout | Includes continuous display of the current interior and high/low temperature, and can be switched from Celsius to Fahrenheit at the touch of a button                                                                                           |
| Large capacity               | Complete 9.88 cu.ft. capacity inside conveniently slim footprint                                                                                                                                                                                |
| High/low temperature alarm   | Audible alarm sounds if the interior temperature rises or falls out of set range                                                                                                                                                                |
| Buffered temperature probe   | Thermometer probe is encased in a glass bead-filled bottle and stored in a fixed box on the center shelf. This is the ideal way to produce temperature readings closest to stored contents, as recommended by the CDC, VFC and VFA's guidelines |
| Access port                  | 3/8" probe hole with cover is pre-installed in unit's side to allow easy installation of user-provided monitoring equipment                                                                                                                     |
| Factory installed lock       | Keyed lock for a secure interior (2 keys included)                                                                                                                                                                                              |
| Automatic defrost            | Reduced user maintenance with automatic defrost                                                                                                                                                                                                 |

## ACR1012GLHD Specifications:

| Overview                   |                     |
|----------------------------|---------------------|
| Nominal Height             | 54.0" (137 cm)      |
| Height of Cabinet          | 53.13" (135 cm)     |
| Nominal Width              | 24.0" (61 cm)       |
| Width                      | 23.38" (59 cm)      |
| Width with Door Open       | 23.88" (61 cm)      |
| Nominal Depth              | 29.0" (74 cm)       |
| Depth                      | 28.33" (72 cm)      |
| Depth with door at 90°     | 50.25" (128 cm)     |
| Capacity                   | 9.88 cu.ft. (280 L) |
| Defrost Type               | Automatic           |
| Door                       | Glass               |
| Cabinet                    | White               |
| US Electrical Safety       | ETL                 |
| Canadian Electrical Safety | ETL-C               |
| Amps                       | 2.0                 |
| Voltage/Frequency          | 115 V AC/60 Hz      |
| Weight                     | 143.3 lbs. (65 kg)  |
| Parts & Labor Warranty     | 1 Year              |
| Compressor Warranty        | 5 Years             |
| UPC                        | 761101107639        |
| Refrigerator Features      |                     |
| Door Swing                 | LHD                 |
| Reversible                 | Factory Reversible  |
| Shelf Type                 | Wire                |
| Shelf Qty                  | 4                   |
| Adjustable Shelves         | Yes                 |
| Thermostat Type            | Digital             |
| Fan Type                   | Interior            |
| Refrigerant Type           | R600a               |
| Refrigerant Amount         | 55.0oz.             |
| High Side PSI              | 0.85                |
| Low Side PSI               | 1.45                |
| Level Legs Qty             | 2                   |
| Interior Light             | Yes                 |
| Temperature Range          | 2° to 8°C           |
| Sound Level                | 37.2 dB             |
| Dimensions                 |                     |

|                        |                 |
|------------------------|-----------------|
| Interior Height        | 44.63" (113 cm) |
| Interior Width         | 20.0" (51 cm)   |
| Interior Depth         | 21.5" (55 cm)   |
| Compressor Step Height | 9.0" (23 cm)    |
| Compressor Step Width  | 20.0" (51 cm)   |
| Compressor Step Depth  | 8.0" (20 cm)    |

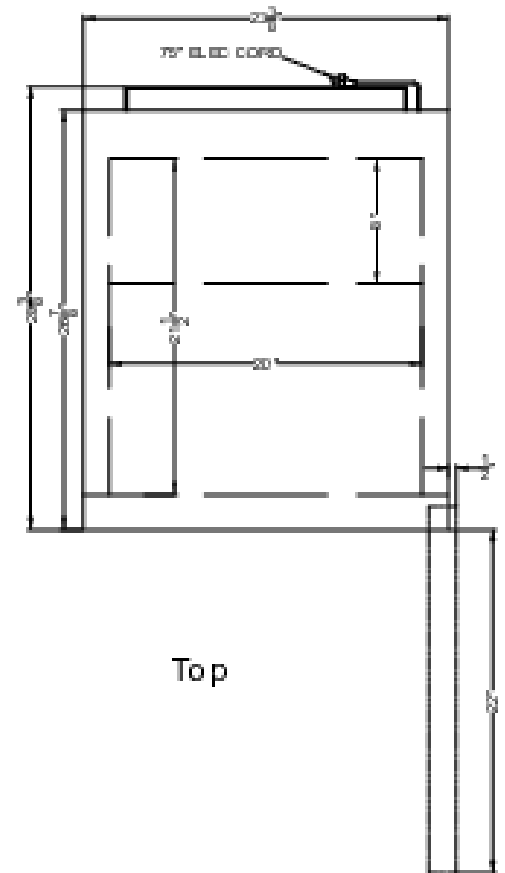

Top

