



 **Convaid**

Convaid Rodeo

Lightweight, Dynamic Positioning, Tilt-in-Space Wheelchair

Convaid Rodeo

The Convaid Rodeo is a lightweight, dynamic positioning, and highly durable wheelchair. The Convaid Rodeo provides up to 45° of adjustable tilt, for individuals that need pressure relief and upper body support. The 90° to 110° recline system opens the seat to back angle for postural relief. The Rodeo is available in four sizes and provides years of growth through its seat depth extension tubes.



Lightweight, Dynamic Positioning and Highly Durable

The Convaid Rodeo wheelchair offers up to 45° adjustable tilt and provides pressure relief and upper body support.



Up to 45° Adjustable Tilt and 110° Recline

Provides pressure relief and upper body support.



Lightweight and Compact-folding

The folding mechanism allows the chair to fold into a tight and lightweight package.



Tension Adjustable Strap Back

Allows for adjustment of the back according to the anatomical shape of the individual to provide optimal support.



Upholstery with Built in Back and Seat Zippers

Allow for foam inserts for additional padding and support.

Standard Features

- 5° to 45° of adjustable tilt
- One-piece height adjustable push handle
- Tension adjustable back
- Cordura upholstery
- Tool-less upholstery removal
- Tool-less adjustable seat depth growth of 5"
- Zippered pockets in seat & back upholstery
- 2-point positioning belt
- Elevating legrest
- One-piece angle adjustable footplates
- 7.5" x 2" front, 11.5" x 2.5" rear solid tires
- Quick release Rolko rear wheels
- Single foot operated wheel lock bar
- Silver powder coated frame
- Closure strap
- Made in the USA

Options

Sizes

(Seat width)

- 10"
- 12"
- 14"
- 16"

Footplates

- One-piece footplate
- Individual footplate

Tires

- 7.5" X 2" front pneumatic tires
- 12.5" rear solid knobby tires

Upholstery Colors

Cordura



Princess pink



Cherry red



Apple green



Forest green



Aqua blue



Royal blue



Sassy purple



Chocolate brown



Panther black

Frame Options



Rodeo with seating



Rodeo without seating (accepts aftermarket seating)

Accessories



Extended headrest cover
- with or without windows
(canopy)



Headrest cover
(canopy)



Rain cover
Canopy required



Mosquito net
Canopy required



Headrest extension



Occi-headwings



Padded headwings



H-harness
(height adjustable)
3-point positioning belt required



Seat reducer



Adjustable firm laterals



Full torso support vest



3-point positioning belt

Accessories



Attendant hand brakes*



Support cushion
(general use)*¹



Position cushion
(medial thigh support/
anti thrust)*¹



Align cushion
(lateral pelvis/lateral thigh
support)*¹



Upper extremity support
surface (clear tray)




Upper extremity support
surface (foam tray)



Medial thigh support
(abductor)



Lateral thigh support
(adductor)

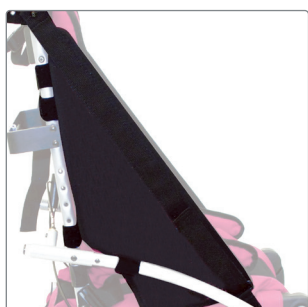


Solid side panels



Mesh side panels



Height adjustable flip-up
armrests

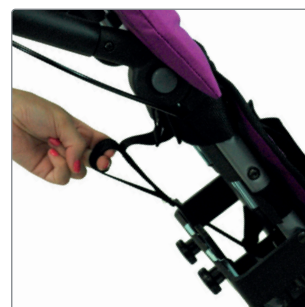

Recline option
(90° to 110°)*

*Manufacturer installation required

¹Incontinence cushion liner available


Accessories



I.V. pole



Oxygen tank bag

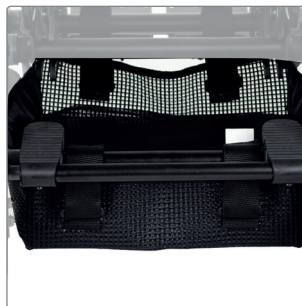


Medical necessity storage basket




Utility bag



Padded elevating legrest cover

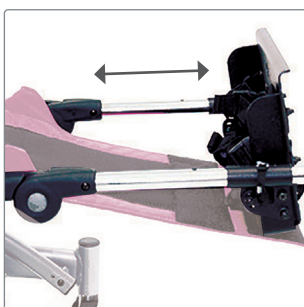

2" elevating legrest extension

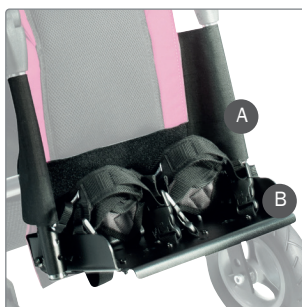

A. Calf panel
B. Foot positioners



Individual angle adjustable footplates



Full padded footbox
One-piece footplate required



Rear anti-tip tubes (pair)



Curb tipper



Castor locks
(for 2" wide front tires)*

*Manufacturer installation required

Specifications

Rodeo measurements		RD 10	RD 12	RD 14	RD 16
	Seat width*	10"	12"	14"	16"
	Seat depth with seating*	7-12"	9-14"	11-16"	13-18"
	Seat depth without seating	9.5-13"	11-15"	13-17"	15-19"
	Seat back height	21"	22"	24"	25"
	Seat to floor	20"	21"	22"	22"
	Recline adjustment	90° - 110°	90° - 110°	90° - 110°	90° - 110°
	Tilt adjustment	5-45°	5-45°	5-45°	5-45°
	Overall height	43"	43"	44"	45"
	Overall length	30.5"	35"	35"	36"
	Overall width	21"	23"	25"	27"
	Folded length	33.5"	34"	34"	40"
	Folded width	21"	23"	25"	27"
	Folded height	17.5"	17"	18"	18"
	Headrest extension height	6"	8"	8"	8"
	Shoulder strap height	13-19"	15-21"	16-22"	18-24"
	Footplate adjustment angle one piece angle adjustable footplate	±20°	±20°	±20°	±20°
	Seat to footplate one piece angle adjustable footplate	9-14" 5.25-11.25**	9-14" 5.25-11.25**	11-16" 7.25-13.25***	13-18" 9.25-15.25**
	Seat to footplate individual angle adjustable footplates	N/A	7-13" 4.5-9.25**	9.5-15.5" 4.25-11***	12-18" 7.25-13.25***
	Push handle adjustment	32-48.5"	34.5-51"	35-52"	36-53.5"
	Weight capacity regular/transit	75/66 lbs	75/66 lbs	100/100 lbs	170/140 lbs
	Weight of chair without front rigging	29 lbs	33 lbs	37 lbs	39 lbs

* Plus or minus .5 inches due to fabric
** Reduced height when footplates are flipped

BR0100-ENG-Rev.00 2022-07-19 | Reserving our rights for printer's error discontinued products. Please find an updated edition on our website © 2022 Etac All rights reserved

For the latest news and continuously updated
product information visit: www.etac.us.com

Etac
+1 844-US-MOBILITY (844-876-6245)
convaidsales.us@etac.com
www.etac.us.com

DCN0952.0


Creating Possibilities