

Delta® Ultra-Light 1000, Full-Electric Bed

HCPCS:



FIGURE A



FIGURE B

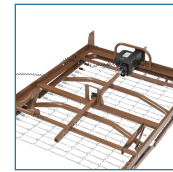


FIGURE C



FIGURE D



FIGURE E

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers' bed ends, old and new
- The transition box mounted on the foot section is key to interchangeability (Figure A)
- The transition box allows you to change the rotation of the high/low shaft
- Headboard is taller than footboard to better fit home decor (Figure B)
- The wood grain panels are break and scratch resistant
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety (Figure D)
- Unique motor is completely self-contained to reduce weight and noise
- The motor assembly can be installed and removed with patient in bed
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility Satisfies all applicable FDA entrapment guidelines
- The motor is UL approved
- In the event of a power failure, a 9 volt battery, located in the hand control, can lower the head and foot sections nine times No crank needed
- New and improved hand pendant has large, easy-to-use controls (Figure E)
- Bed ships in two cartons
- Channel frame construction provides superior strength and reduced weight
- Zinc-coated spring deck
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84"
- Patent Pending
- Weight Capacity 450 lbs

PRODUCTS	Description
15033	UltraLightFullElec.Bed
15033BV-FR	UltraLightFullElec.Bed w/ FullRails
15033BV-HR	UltraLightFul Elec.Bed w/ HalfRails
15033BV-PKG	UltraLightFullElec.Bed FR Pkg Spr. Matt
15033BV-PKG-1	UltraLightFullElec.Bed HR Pkg Spr. Matt
15033BV-PKG-11	UltraLightFullElec.Bed FR Pkg QNE Matt
15033BV-PKG-12	UltraLightFullElec.Bed HR Pkg QNE Matt
15033BV-PKG-1-T	UltraLightFullElec.Bed HR Pkg TSS Matt
15033BV-PKG-2	UltraLightFullElec.Bed FR Pkg Foam Matt
15033BV-PKG-4	Full Elect. Ultra Lt. Plus Bed,1/cs
15033BV-PKG-T	UltraLightFullElec.Bed FR Pkg TSS Matt

SPECIFICATIONS

	Description
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components
Weight Capacity	450 lbs Maximum weight limit includes weight of occupant
Carton Shipping Weight	182 lbs.
Weight	169 lbs.
CA Prop 65	[WARNING: This product can expose you to chemicals that are known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov .]

ACCESSORIES

	Description
15030HBL	Hospital Bed Linen Kit,1 EA,10/CS, CS
15030HBC	Bedding in a Box,1 EA,5/CS, CS