

Gravity 9 Premium Long Term Care Pressure Redistribution Mattress

HCPCS:

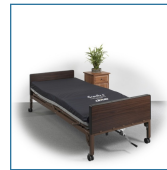


FIGURE A

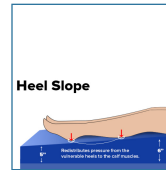


FIGURE B

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- Drive's premium therapeutic Gravity foam mattress maximizes value and comfort with multiple layers of latex free foam
- Full top layer of open cell, visco-elastic memory foam improves immersion and pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with channel cut foam center offloading pressure from vulnerable heels
- Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

PRODUCTS	Description
15970	Gravity 9 Mattress 80", 1/EA
15970-4280	Gravity 9 Mattress, 42"W x 80"L, 1 EA
15970-4280-RR	Gravity 9 Matt , 42"W x 80"L w/RR, 1
15970-4284	Gravity 9 Mattress, 42"W x 84"L, 1 EA
15977	Gravity 9 MattressW/Elev Perim 80", 1/EA
15984	Gravity 9 Mattress 84", 1/EA
15996	Gravity 9 Mattress 76", 1/EA
G9-39X80X6	Gravity 9 Mattress 39"W x 80"L x 6"H

SPECIFICATIONS	Description
Warranty	36" W -15 Years Limited Foam 3 Year Cover 42" W- 7 Years Limited Foam 3 Year Cover
Weight Capacity	450 lbs
Side Rails	Yes
Cover Material	Stretch, Fluid-Resistant Nylon Top Waterproof Vinyl Bottom
CA Prop 65	[WARNING: This product can expose you to chemicals that are known to the State of California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov .]

ACCESSORIES	Description
15977C	Mattress Cover 15977C,1/ea, EA
15970C	Mattress Cover 15970C,1/ea, EA
15984C	Mattress Cover 15984C,1/ea, EA
15996C	Mattress Cover 15996C,1/ea, EA