

Balanced Aire Non-Powered Self Adjusting Convertible Mattress

HCPCS: E0373, E0277,E0373, E0277



FIGURE A

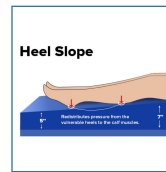


FIGURE B



FIGURE C

- Versatile reactive mattress replacement system with self adjusting air/foam cells provides outstanding pressure redistribution and shear/friction reduction
- Designed to prevent, treat and heal pressure ulcers in the home or long term care setting
- Add an alternating pressure pump and the Balanced Aire easily converts to a dynamic powered system (14005E or BA9600-PUMP)
- 2.5" Hybrid Visco foam top layer provides immersion, premium comfort and heat dissipation
- High density foam perimeter provides firm edge support and facilitates safe patient transfers
- Heel slope redistributes pressure away from vulnerable heels
- Self-adjustment can be open cell or closed cell to meet individual needs
- Self-Adjusting air cells recognize body weight and pressure and adjust to reduce interface pressure and improve patient comfort
- Eight non-circular cylinders prevent "rolling" or "gapping"
- Removable 4 way stretch PU cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Arrives patient ready easy to deliver, set up and transport
- Foam and cover are antimicrobial, antifungal, and latex free
- Meets flammability standards CFR1632 and CRF1633

PRODUCTS	Description
BA9600-NP	35X80X7 Non-Pwr SelfAdjAir/Foam Matt,1/e
BA9600-NP-42	42X80X7 Non-Pwr, SelfAdjAir/FoamMatt,1/e
BA9600-NP-84	35X84X7 NonPwr SelfAdjAir/Foam Matt, 1/e
BA9600-P	36X80X7 PwrdSELF ADJ.AIR/FOAM,MATT,1/E
BA9600-P-42	42X80X7 Pwrd Self Adj.Air/Foam,Matt,1/e
BA9600-P-84	36X84X7 PwrdSelfAdjAir/Foam Matt, 1/e

SPECIFICATIONS	Description
Warranty	5 Year Limited Mattress 1 Year Control Unit
Weight Capacity	35" W 500 lbs
Control Unit Airflow	When Used w/Power Source (P) 4 LPM
Control Unit Alarms	N/A
Control Unit Cycle Time	When Used w/Power Source (P) 5 Minutes
Control Unit Type	When Used w/Power Source (P) Compressor
Mattress Cover Material	Fluid-Resistant Stretch

ACCESSORIES	Description
BA9600-C-NP3580	Replacement Cover, BA9600-NP(P),1 ea, EA, EA

ACCESSORIES**Description**

BA9600-P-42-TC	Replacement Cover,BA9600-NP(P)-42,1ea,EA, EA
BA9600-C-NP3584	Replacement Cover,BA9600-NP(P)-84,1ea,EA, EA
14005E	Pump only 14001E, 1/EA, EA
14427QF	Quick Fill Pump w/adapters, 1 EA/12 CV, EA