

Delta® Ultra-Light 1000, Full-Electric Bed

HCPCS:



FIGURE A



FIGURE B

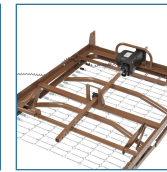


FIGURE C



FIGURE D



FIGURE E

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers' bed ends, old and new
- The transition box mounted on the foot section is key to interchangeability (Figure A)
- The transition box allows you to change the rotation of the high/low shaft
- Headboard is taller than footboard to better fit home decor (Figure B)
- The wood grain panels are break and scratch resistant
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety (Figure D)
- Unique motor is completely self-contained to reduce weight and noise
- The motor assembly can be installed and removed with patient in bed
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility Satisfies all applicable FDA entrapment guidelines
- The motor is UL approved
- In the event of a power failure, a 9 volt battery, located in the hand control, can lower the head and foot sections nine times No crank needed
- New and improved hand pendant has large, easy-to-use controls (Figure E)
- Bed ships in two cartons
- Channel frame construction provides superior strength and reduced weight
- Zinc-coated spring deck
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84"
- Patent Pending
- Weight Capacity 450 lbs

PRODUCTS	Description
15033	Full Electric Ultra Light Plus,1/cs
15033BV-FR	Ultra Lt.Full Elec.Bedw/Full R
15033BV-HR	Ultra Lt.Full Elec.Bed,w/Half
15033BV-PKG	Ultra Lt.Full Elec.Bed Pkg,1pk
15033BV-PKG-1	Ultra Lt.Full Bed Pkg w/HalfRa
15033BV-PKG-11	Delta Full,FullRail,QuickEasyFoam,1ea
15033BV-PKG-12	Delta Full,HalfRail,QuickEasyFoam,1ea
15033BV-PKG-1-T	Ultra Lt.Full Bed Pkg w/BrownV
15033BV-PKG-2	Ultra Lt.Full Elec.BedPkg w/Fo
15033BV-PKG-4	Full Elect. Ultra Lt. Plus Bed,1/cs
15033BV-PKG-T	Ultra Lt.Full Elec.Bed Pkg,1pk

SPECIFICATIONS**Description**

Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components
Weight Capacity	450 lbs Maximum weight limit includes weight of occupant
Carton Shipping Weight	182 lbs.
Weight	169 lbs.

ACCESSORIES**Description**

15030HBL	Hospital Bed Linen Kit,1 EA,10/CS, CS
15030HBC	Bedding in a Box,1 EA,5/CS, CS