



LED controls with optional keypad

Model configurations			
Controls	Keypad	V/Hz/Ph	Item number
LED	No	115/60/1	<b>FZR1</b>
	Yes	115/60/1	<b>FZR1-KP</b>

## Options and accessories<sup>1</sup>

- Replacement NTC probe with NIST-traceable certificate of calibration (item# 01077841)
- Wall mount bracket (item# 00969949)
- Organizing bins, 6 (item# 00969931)
- Extra shelf (item# 00969956)
- Temperature surveillance module (includes alarm, 8.00" (20.32 cm) chart recorder, remote contacts) (item# 01570860)
- Glycerine, 16 oz (item# 00959296)
- Casters (set of 4) (item# 01111947)
- Pedestal, 2.00" (5.08 cm) (item# 01184845)
- Pedestal, fan-assist, 4.50" (11.43 cm) (item# 01025899)
- Pedestal, 12.00" (30.5 cm) (item# 01025881)
- NO/NC dry contacts for connection to remote alarm systems (item# 01092022)
- Stacking kit (countertop on countertop) (item# 01100742)
- Stacking kit (countertop on Performance Plus undercounter) (item# 01100759)
- Stacking kit (countertop on REF/FZR (Legacy) undercounter) (item# 01142363)

<sup>1</sup> Certifications do not apply to options and accessories.

**FZR1**

medical-grade countertop freezer

## Features

### Superior temperature performance

- heavy-duty compressor and forced air cooling
- quiet compressor and coaxial fans
- non-CFC R404a refrigerant
- automatic defrost

### Convenience and security

- exterior LED displays product temperature in user-selectable F or C
- "sleep mode" turns off display
- user-programmable, audible and visual high/low alarms
- lockable controller prevents inadvertent setpoint changes
- adjustable refrigeration set points are accessed through front controls
- high/low temperature log shows highest and lowest product temperature recorded
- cylinder door lock
- back wall port provides entry for third-party probes
- flush door design is compatible with standard automated dispensing system brackets

### Durability and serviceability

- stainless steel construction on exterior and interior
- rubberized feet
- epoxy-coated wire shelf
- heavy-duty edge-mount self-closing hinges
- field-reversible door
- removable back panel for easy service access to compressor, condenser and fans
- dart-type magnetic gasket
- dished floor

### Warranty

- 2 year parts and labor
- 5 year compressor parts only
- optional extended one year warranty (item# EW12)
  - 3 year parts and labor
  - 5 year compressor parts only
- optional extended three year warranty (item# EW36)
  - 5 year parts and labor

### Certifications (applied to all products unless otherwise noted)



Commercial Refrigerator  
SAFETY US/CA – SA12646

## Specification

Nominal capacity	0.6 cu ft (16 L)
Ventilation clearance with more than 12.00" (30.5 cm) top clearance	2.00" (5.1 cm) for each side
Ventilation clearance with less than 12.00" (30.5 cm) top clearance	6.00" (15.3 cm) for one side

## Exterior

W1 Width	18.60" (47.2 cm)
Depth at top	24.00" (61.0 cm)
D1 Depth at bottom	20.00" (50.8 cm)
H1 Height	17.93" (45.5 cm)

Door swing	standard with self-closing right hinge, field reversible
------------	--

## Interior

Width	14.60" (37.1 cm)
Depth	11.60" (29.5 cm)
Height	5.72" (14.5 cm)
Storage system	(1) 13.81" x 9.06" (35.08 cm x 23.01 cm), epoxy-coated wire shelf

## Door configurations

Door	stainless steel
Door lock	cylinder door lock with 2 keys
Door hinges	edge-mount, self-closing
Gasket	magnetic, dart-style

## Electrical

C1 115 V/60/1 electrical	4.1 full load amps, NEMA 5-15P 90° hospital-grade plug low profile, SJT, 18 AWG, 8' (2.4 m) cord
--------------------------	--

## Controls

Temperature display	standard 7 segment LED, default set in C, user-selectable F or C
Display sleep function	button turns off display, press again to turn back on
Factory preset	-29 C
Product simulation bottle	(1) 60 ml
Probes	3 probes : P1 refrigeration system; P2 defrost; P3 bottle
Programmable operating range	-28 to -18 C (-20 to 0 F)

## Alarming

High and low product temperature alarm	audible and visual, user programmable set points in F or C
Probe error alarm	audible and visual — product temperature probes, refrigeration and defrost
Alarm delay	user-programmable 0 to 60 minute (default 1 minute)
Alarm mute	user selectable
Alarm volume	0-10 adjustable

## System performance

Refrigerant and charge	R404a refrigerant, 9.2 oz
Automatic defrost	electric, timed, adjustable
Energy consumption	5.4 kWh/day
Max heat rejection	875 BTU/hr

## Shipping



Approximate net weight	72 lb (33 kg)
Approximate ship weight	108 lb (49 kg)

NOTE: For indoor use only

## FZR1 countertop freezer

**SHORT FORM SPECIFICATION:** Freezer with solid stainless steel door fits on standard 24.00" (61.0 cm) deep countertop and under 18.00" (45.7 cm) wall hung cabinets with 0.6 cu ft (16 L) of nominal capacity. Includes (1) epoxy-coated wire shelf. Environmentally responsible R404a forced-air cooled refrigeration system. Bottom-mounted controls display product temperature in user selectable C or F. Integral high and low temperature alarming. Storage area insulated with CARB compliant non-HFC foam. 8' (2.4 m) power cord with NEMA 5-15p 90 degree hospital-grade plug. UL and cUL listed.

## Dimensional drawing

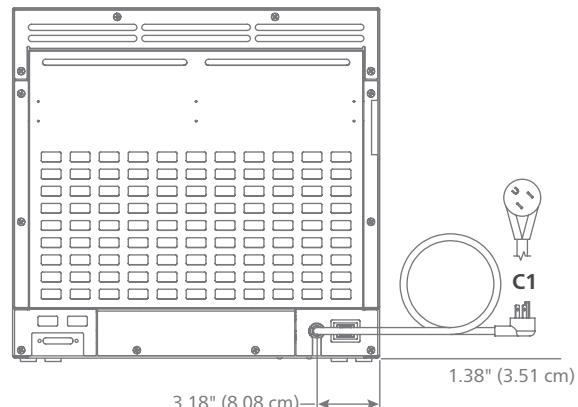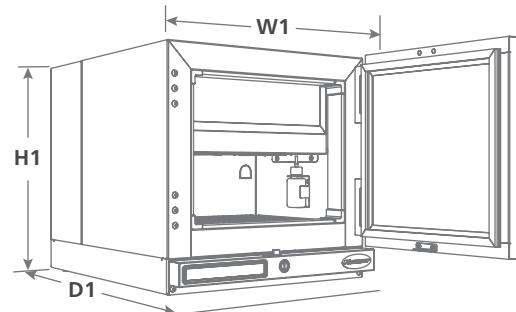

FOLLETT and FOLLETT HEALTHCARE are registered trademarks of Follett Products, LLC, registered in the US.

Follett reserves the right to change specifications at any time without obligation.

Certifications may vary depending on country of origin.