


SFPvh 1402

50 cu.ft. Perfection High performance Laboratory Freezer

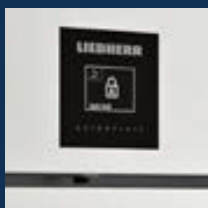


SmartLock

To prevent products from accidentally entering circulation after an interruption in the cold chain, SmartLock automatically locks the door if there is an alarm. Only authorized persons can unlock it again.




Integrated WiFi/LAN interface

The integrated WiFi/LAN interface enables you to easily connect your appliances with existing local systems.

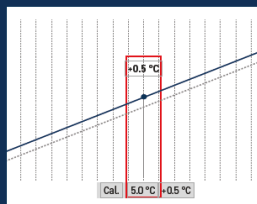

Calibration function

Due to operational requirements and application-specific conditions, it is possible to set temperature deviation of the interior from the set temperature to a storage temperature determined with calibrated thermometers. These deviations can be adjusted by the calibration function. For the secure knowledge of cooling at precise temperatures.

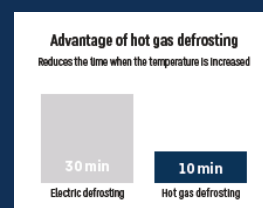

Energy-efficient hot gas defrosting


Timed and on-demand hot gas defrosting allows for very short defrosting times. To further improve temperature stability during defrosting, the internal temperature is slightly lowered shortly beforehand to create a reserve of cold air.

LIEBHERR

SFPvh 1402

50 cu.ft. Perfection High performance
Laboratory Freezer

Laboratory Freezer

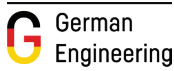

home.liebherr.com

SFPvh 1402


Total gross/net capacity l (cu.ft.)	1391 / 1091 l (48.2 / 46.8 cu. ft.)
Exterior dimensions in mm (w/d/h) in inches (w/d/h)	1427 / 867 / 2135 mm 56 1/4" / 34 3/16" / 84 1/16"
Interior dimensions in mm (w/d/h) in inches (w/d/h)	1210/638/1493 mm 47 11/16" / 25 1/8" / 58 13/16"
Energy Consumption (kWh/day) *	5.26 kWh
EnergyStar rated	Yes
Ambient temperature in °C (°F)	+10°C to +35°C (50°F to 95°F)
Refrigerant	R 290
Sound rating**	60 dB(A)
Connection rating	4.0 A
Voltage/ Frequency	115 V / 60Hz
Cooling system/Defrost	Forced air / automatic
Temperature range	-9°C to -35°C
External cabinet finish	Steel/ white
Door/lid material	Steel/ white
Internal liner material	Chrome-nickle steel
Type of control	2.4" touch & swipe color display
Temperature display	External, digital
Alarm	Visual, audible
Alarm forwarding	Yes
Connectivity solution / Interface	Integrated, removable / WiFi/LAN
Data Logger	integrated, USB interface
Interior Light	LED
Adjustable shelves	8
Shelf material	Wire shelves plastic-coated
Shelf load weight in kg freezer	60 kg
Usable Shelf Area in mm (inches)	480/638 (24 1/16"/19 11/16")
Handle	Vertical recessed handle
Lock	Electronic
Self-closing door	Yes
Door Hinges	Left- hand/ right hand
Gross/net weight in Kg (lbs)	235/210 (518/463)

* Measured unloaded at 26.6°C (80°F) ambient temperature with set point of 3.3°C (38°F).

** Noise level – sound power re 1pW – according to IEC 60704-2-14.



Unit Dimensions

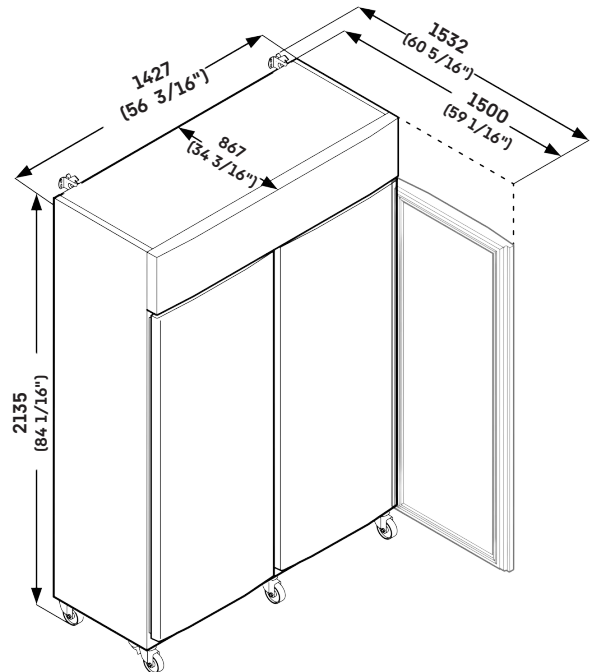

Accessories

U-shaped trayslide, (left / right)	7745088 / 7745090
Gray wire shelf	7112608
Castors Ø 100 mm (4 in)	9597965
Adjustable feet, 115 – 175 mm (4.5 - 6.9 in)	9080132
Foot pedal opener	9597649
Side-by-Side set in white	9597677

LIEBHERR