


SFTsg 1501

3.6 cu. ft. Performance undercounter freezer

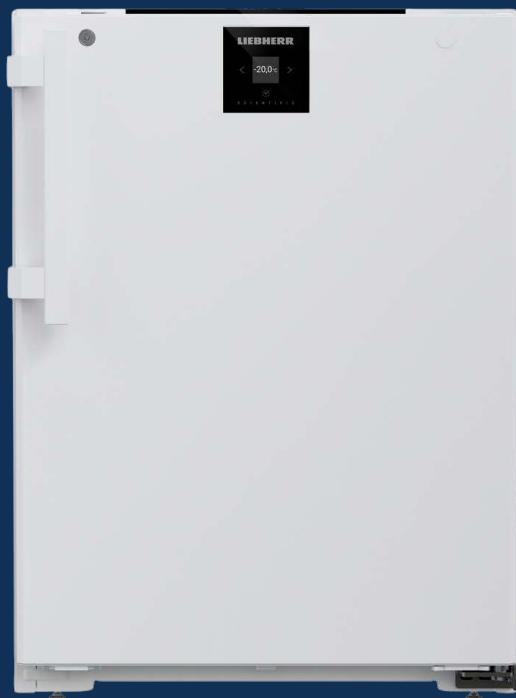


Energy Star

The LED front lighting provides optimal illumination of the drawers even in dark environments. Thanks to the perfect overview, the desired substance is immediately at hand, the door can be quickly closed again, and the interior temperature remains stable. For low consumption and lower energy costs.




Retrofittable WiFi/LAN interface

Energy-efficient refrigerators and freezers help with ever increasing electricity prices. Liebherr fully utilizes the savings potential of every component – from the compressor to the seal, to the lighting and electronics. The ultra-efficient appliances in the Perfection series for example demonstrates how this has a positive effect on your overall operating costs.

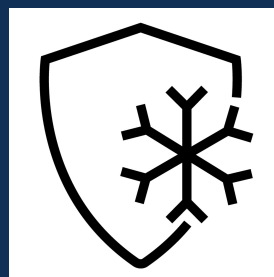

SmartFrost

SmartFrost reduces frost build-up in the interior and on stored products. Defrosting is required less often – and when it is necessary, it is made much easier by the smooth, easy-to-clean interior walls. The evaporator in the foam insulation ensures high levels of energy efficiency as well as even cooling.

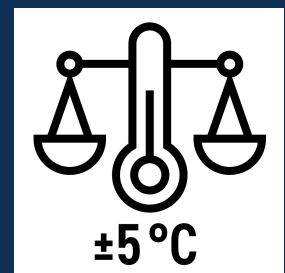

±5 °C temperature stability


Sensitive substances cannot be subjected to large deviations from the set storage temperature. Liebherr freezers guarantee a temperature fluctuation of no more than ±5 °C, so that the valuable goods you store in them are optimally protected and their quality is preserved.

LIEBHERR

SFTsg 1501

3.6 cu. ft. Performance
undercounter freezer


Laboratory Freezer

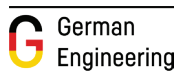

home.liebherr.com



Unit Dimensions

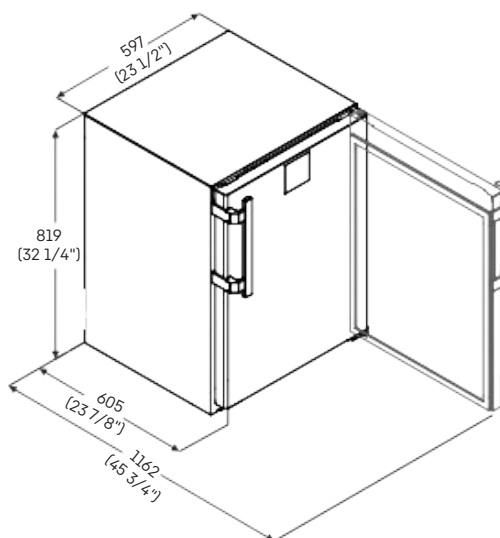

Accessories

Set of 3 drawers	987626
Castors / Adjustable legs	9904819 / 9906105
Additional lock barrels (up to 5)	on request
Stacking kit, white	9876098
Side-by-side set in white	9031901
SmartModule	6145269

SFTsg 1501	
Total gross/net capacity l (cu.ft.)	3.6 / 3.5 cu. ft. (101 / 99 l)
Exterior dimensions in mm (w/d/h) in inches (w/d/h)	597 / 605 / 819 23 1/2" / 23 7/8" / 32 1/4"
Interior dimensions in mm (w/d/h) in inches (w/d/h)	380 / 410 / 610 15" / 16 1/8" / 24"
Energy Consumption*	0.49 kWh / day
EnergyStar rated	Yes
Ambient temperature in °C (°F)	+10°C to +35°C (+50°F to +95°F)
Refrigerant	R 600a
Sound rating**	40 dB(A)
Connection rating	2.0 A
Voltage / Frequency	115V / 60 Hz
Cooling system/Defrost	SmartFrost / manual
Temperature range	-9 °C to -30 °C
External cabinet finish	Steel / white
Door/lid material	Steel / white
Internal liner material	Commercial grade moulded polystyrene
Type of control	Monochrome display, touch buttons
Alarm	Visual and audible
Alarm forwarding	Yes
Interior Light	-
Adjustable shelves	3
Shelf material	wire shelves plastic-coated
Shelf load weight	45 kg (99 lbs)
Usable Shelf Area	380x 338 mm (15" x 13 1/4")
Handle	Antimicrobial handle with opening mechanism
Lock	Mechanical
Self-closing door	Yes
Door Hinges	Right, reversible on site
Access port	1 x Ø 10 mm
Gross/net weight	37 / 34 kg (81 / 75 lbs)

* Measured unloaded at 26.6°C (80°F) ambient temperature with set point of 3.3°C (38°F).

** Noise level – sound power re 1pW – according to IEC 60704-2-14.

V.O (June 2025)

LIEBHERR