


# SRPvh 6511

21 cu. Ft. Perfection High-performance Laboratory Refrigerator with glass door



## Integrated WiFi/LAN interface

The integrated WiFi/LAN interface enables you to easily connect your appliances with existing local systems.

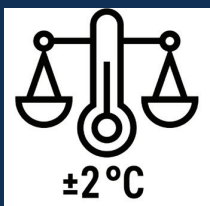

## $\pm 2^{\circ}\text{C}$ temperature stability

Temperature fluctuations must be minimal to ensure that the quality of sensitive substances is preserved. Liebherr appliances have an impressive maximum deviation of just  $\pm 2^{\circ}\text{C}$  from the set temperature. This not only maintains the storage climate, but also the value of your stored goods.



## SmartLock

To prevent products from accidentally entering circulation after an interruption in the cold chain, SmartLock automatically locks the door if there is an alarm. Only authorized persons can unlock it again.



## Vertical LED lighting


Vertical LED lighting provides optimal illumination of the interior over the entire height of the appliance. This not only saves you time looking for things, but, above all, cuts energy costs.

# LIEBHERR

# SRPvh 6511

21 cu. Ft. Perfection High-performance  
Laboratory Refrigerator with glass door

## Laboratory Refrigerator

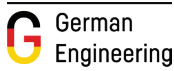

home.liebherr.com

### SRPvh 6511

Total gross/net capacity	642 / 480 l (22.2 / 20.7 cu. ft.)
Exterior dimensions in mm (w/d/h) in inches (w/d/h)	697 / 867 / 2135 27 1/2" / 34 3/16" / 84 1/16"
Interior dimensions in mm (w/d/h) in inches (w/d/h)	489 / 638 / 1493 19 5/16" / 25 1/8" / 58 13/16"
Energy Consumption*	1.01 kWh/day
EnergyStar rated	Yes
Ambient temperature	+10°C to +35°C (50°F to 95°F)
Refrigerant	R 600a
Sound rating**	60 dB(A)
Connection rating	2.0 A
Voltage / Frequency	115 V / 60Hz
Cooling system/Defrost	Dynamic
Temperature range	0°C to +16°C
External cabinet finish	Steel/white
Door/lid material	Insulating glass door with steel frame
Internal liner material	Chrome-nickle steel
Type of control	2.4" touch and swipe color control
Temperature display	External digital
Connectivity Solution/ Interface	Integrated, removable/WIFI/LAN
Data logger	Integrated, USB interface
Alarm	Visual, audible
Alarm forwarding	Yes
Interior Light	Vertical LED light strip
Adjustable shelves	4
Shelf material	Wire shelves plastic-coated
Shelf load weight	60 kg (132 lbs)
Usable Shelf Area	480/638 mm (18 7/8" / 25")
Handle	antimicrobial handle with opening mechanism
Lock	Electronic
Self-closing door	Yes
Door Hinges	Right, reversible on site
Insulation	70 mm (2 13/16")
Gross/net weight	220/201 kg (489.4/443.1 lbs)


\* Measured unloaded at 26.6°C (80°F) ambient temperature with set point of 3.3°C (38°F).

\*\* Noise level – sound power re 1pW – according to IEC 60704-2-14.

V.0 (May 2025)



### Unit Dimensions

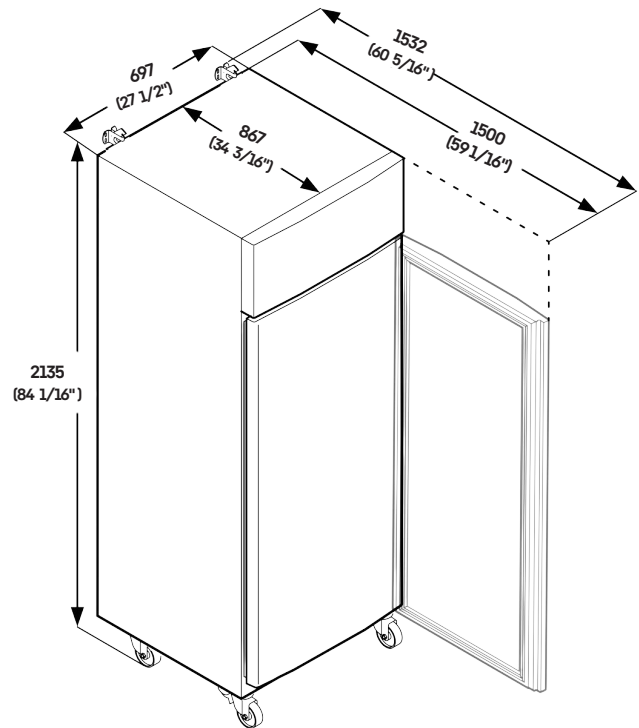

### Accessories

U-shaped trayslide, (left / right)	7745088/ 7745090
Gray wire shelf	7112608
Castors Ø 100 mm (4 in)	9597828
Adjustable feet, 115 – 175 mm (4.5 - 6.9 in)	9080131
Foot pedal opener	9597649
Side by side set in white	9597677

# LIEBHERR