
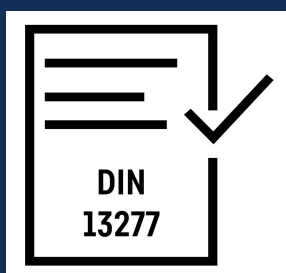


SRTvh 1501

4.7 cu. ft. Perfection solid door undercounter refrigerator

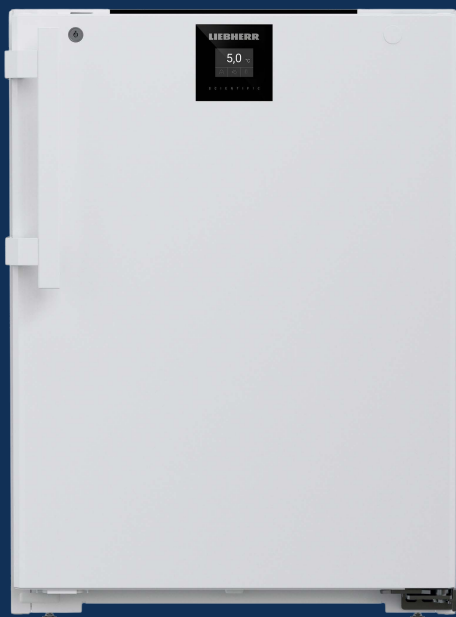

DIN 1277 compliant

Laboratory and pharmacy refrigerators from Liebherr meet all the necessary features and performance criteria to fully comply with DIN 13277. DIN 13277 defines the requirements and test specifications for refrigerators and freezers used in laboratories and for storing medicinal products.




Ultra-efficient

Energy-efficient refrigerators and freezers help with ever increasing electricity prices. Liebherr fully utilises the savings potential of every component – from the compressor to the seal, to the lighting and electronics. The ultra-efficient appliances in the Perfection series for example demonstrates how this has a positive effect on your overall operating costs.



+2 °C SafetyDevice

If the temperature falls below +2 °C due to a malfunction while the temperature setting is at +5 °C or higher, the SafetyDevice takes control, stabilizing the temperature, activating the alarm and protecting the contents from freezing. At colder settings, this activation threshold decreases relatively with the temperature setting.

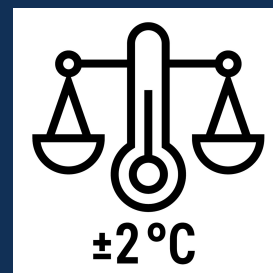

±2 °C temperature stability


Temperature fluctuations must be minimal to ensure that the quality of sensitive substances is preserved. Liebherr appliances have an impressive maximum deviation of just ±2 °C from the set temperature. This not only maintains the storage climate, but also the value of your stored goods.

LIEBHERR

SRTvh 1501

4.7 cu. ft. Perfection solid door undercounter refrigerator


Laboratory Refrigerator

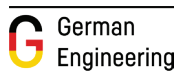

home.liebherr.com



Unit Dimensions

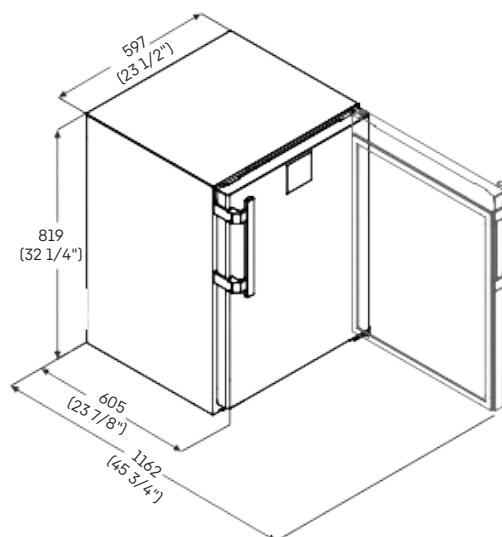

SRTvh 1501	
Total gross/net capacity	4.7 / 4.0 cu. ft. (133 / 113 l)
Exterior dimensions in mm (w/d/h) in inches (w/d/h)	597 / 605 / 819 23 1/2" / 23 7/8" / 32 1/4"
Interior dimensions in mm (w/d/h) in inches (w/d/h)	450 / 440 / 620 17 3/4" / 17 3/8" / 24 3/8"
Energy Consumption*	0.34 kWh/day
EnergyStar rated	Yes
Ambient temperature	+10°C to +35°C (+50°F to +95°F)
Refrigerant	R 600a
Sound rating**	43 dB(A)
Connection rating	2.0 A
Voltage/ Frequency	115V / 60 Hz
Cooling system/Defrost	Forced-air / automatic
Temperature range	+3°C to +16°C
External cabinet finish	Steel / white
Door/lid material	Steel / white
Internal liner material	commercial grade moulded polystyrene
Type of control	Touch & Swipe
Alarm	Visual and audible
Alarm forwarding	Yes
Interior Light	LED lighting column
Adjustable shelves	3
Shelf material	Wire shelves plastic coated
Shelf load weight	45 kg (99 lbs)
Usable Shelf Area	450 x 425 mm (17 3/4" x 16 3/4")
Handle	Antimicrobial handle with opening mechanism
Lock	Mechanical
Self-closing door	Yes
Door Hinges	Right, reversible
Access ports	1 x Ø 12 mm
Gross/net weight	37 / 33 kg (81.5 / 72.75 lbs)

Accessories

White grid shelf	7113408
Perforated metal shelf, white	7746588
Castors/ adjustable legs	9904819/ 9906105
Scanner rail	7401607
Additional lock barrels	on request
Stacking kit, white	9876098
Side-by-side, set in white	9031901
SmartModule	6145269

* Measured unloaded at 26.6°C (80°F) ambient temperature with set point of 3.3°C (38°F).

** Noise level – sound power re 1pW – according to IEC 60704-2-14.

V.0 (June, 2025)

LIEBHERR