



24" ADA Height General Purpose All-Refrigerator

Model # MSRE224



Standard Features

- 1" diameter probe port in rear cabinet wall (standard)
- 3 full-depth perforated metal shelves, cantilever mounted on rear wall of interior, adjustable in 1.5" increments
- Flat interior door liner maximizes depth of shelves
- Front vented through toe space grill for true built-in installations
- Extra-long, jacketed commercial power cord with hospital grade plug to meet stringent hospital and clinical use standards
- UL Approved and warrantied for commercial use
- Standard door lock for secure storage needs
- Discrete top edge mounted door handle
- Cam-assisted door close hinge system to reduce door ajar events
- Bright energy efficient white LED interior lighting for easy content viewing
- BrightShield® Antimicrobial Lighting Technology creates a more hygienic environment
- 4 leveling legs for up to 1" height adjustment
- Durable white powder coat paint finish
- Recommended ambient temperature 18.3°C to 32.2°C (65°F to 90°F)
- Additional shelves can be ordered separately, the ordering part number is S41013646-WHT

Reliable Temperature Stability & Energy Efficiency

- Advanced forced air evaporator system design with optimized interior cooling air distribution
- Powerful compressor for fast pull-down time and temperature recovery after door opening or content loading
- Automatic defrost without supplemental heat

Temperature Control

- Door-front mounted control for easy access and monitoring
- Accurate electronic temperature control of contents within $\pm 2^{\circ}\text{C}$ of set point
- Large, easy-to-read digital temperature display
- User-selectable temperature set point between $+2^{\circ}$ and $+8^{\circ}\text{C}$, factory preset to $+4^{\circ}\text{C}$
- User-selectable high and low temperature alarm values
- Visual and audible alarms for over temperature, under temperature, and door ajar events
- Visual diagnostic error codes

Rough-In Dimensions		
"A"	"B"	"C"
24" (61 cm)	31 1/4" to 32 1/4" 79.4 to 81.9 cm)	*

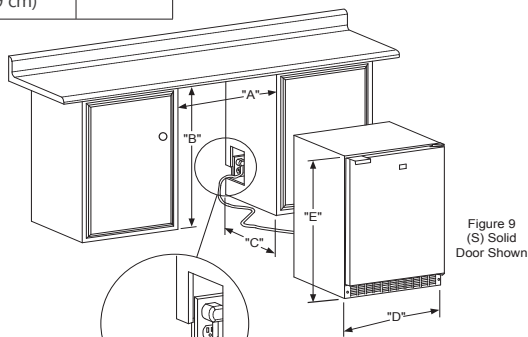


Figure 9a



Figure 10

If necessary to gain clearance inside the rough-in opening a hole can be cut through the adjacent cabinet and the power cord routed through this hole to a power outlet. Another way to increase the available opening depth is to recess the power outlet into the rear wall to gain the thickness of the power cord plug. Not all recessed outlet boxes will work for this application as they are too narrow, but a recessed outlet box equivalent to Arlington #DVR1W is recommended for this application, (see Figure 10).

Storage Configuration

(3) cantilever, fully adjustable white steel shelves

Product Dimensions

23 7/8" W x 31 1/4" minimum H x 23 23/32" D
(excluding door handle)

Finish/Door Options




White



White Framed Glass

Cabinet Dimensions					
"D"	"E"	"F"	"G"	"H"	"J"
23 7/8" (60.7 cm)	31" to 32" 79 to 81 cm	23 23/32" (60.2 cm)	24 3/32" (63.4 cm)	46 3/8" (118 cm)	24 7/8" (63.2 cm)

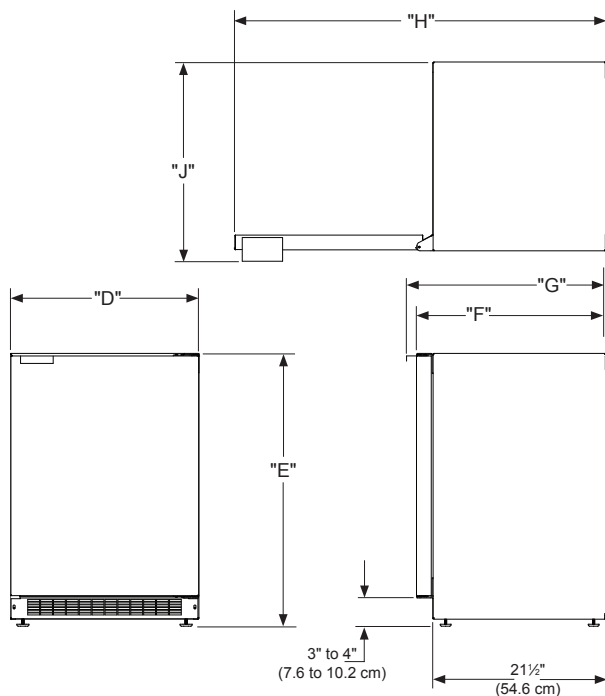

Ordering Details

Solid White Door	MSRE224-WS81A
White Framed Glass Door	MSRE224-WG81A

Width	23 7/8"
Interior Finish	White
Capacity	4.9 Cubic Feet 139 Liters
Shelving System	Cantilever
Shelf Fronts	White
Interior Lighting	White LED + BrightShield
Toe Grill Finish	Black
Hinge	Cam-assisted self close

Handle	Discrete
Control Type	Electronic
Lock	Yes
Temperature Range	36°F to 46°F (2°C to 8°C)
Electrical Requirements	120V / 60Hz / 15A
Length of Power Cord	5'

Product Dimensions	23 7/8" W x 31" H x 23 23/32" D
Depth to Front of Door	23 23/32"
Depth to Handle	24 31/32"
Depth with door at 90°	46 3/8"
Shipping Weight	135 lbs
Finishes	Solid white, white framed glass
Agency Approvals	UL