

Silicone Adhesive Foam Dressings

McKesson silicone adhesive foam dressings are absorbent and atraumatic dressings that absorb exudate, maintain a moist wound-healing environment, and act as a barrier to liquids and microorganisms from the wound. The wound contact surface of the silicone adhesive foam dressing is coated with a layer of soft silicone that does not adhere to the wound surface or cause trauma to delicate new tissue upon removal. This soft silicone layer is also slightly tacky, which facilitates application and retention of the dressing to intact skin, but does not cause epidermal stripping or pain on removal.

Indications: McKesson silicone adhesive foam dressings are indicated for the management of exuding wounds, such as leg and foot ulcers, pressure ulcers, traumatic and surgical wounds, as well as superficial and partial-thickness burns. The dressing absorbs exudate and maintains a moist wound-healing environment while minimizing the risk of maceration. It may be applied and changed by clinicians or patients.

Contraindications: Discontinue use if the patient is allergic to any components of the dressing.

Change Frequency: The interval between dressing changes will depend entirely upon the state of the wound and the nature of the secondary dressing. Dressing applied to heavily exuding or sloughy wounds may need to be replaced daily initially; however, as healing progresses and the amount of exudate decreases, the interval between changes can be extended to up to 7 days or as directed by healthcare professional.

Features and Benefits

- Perforated silicone gel adhesive provides gentle but firm adhesion
- Highly conformable and remains intact on skin
- Soft silicone minimizes pain to patients and trauma to wounds upon removal
- Affixes securely to dry, intact skin; seals the wound margins to help prevent exudate from leaking on surrounding tissue, minimizing risk of maceration
- Breathable, moisture and bacteria resistant outer film layer
- Absorbent inner hydrocellular foam pad maintains a moist wound environment to help promote healing
- Dressing may be lifted and repositioned without losing its adherent properties
- Wear time up to 7 days or as directed by a healthcare professional
- Not made with natural rubber latex




Silicone Adhesive Foam Dressings



Description	E1 #	MFR #	QTY	HCPCS
Hydrocellular Foam Dressing, With Border, Silicone Adhesive, Sterile				
3 in x 3 in, 5-Layer	886431	4842	1/ea, 10/bx, 20 bx/cs	A6212
4 in x 4 in, 5-Layer	886432	4843	1/ea, 10/bx, 20 bx/cs	A6212
5 in x 5 in, 5-Layer	1276384	16-4664	1/ea, 10/bx, 20 bx/cs	
6 in x 6 in, 5-Layer	886433	4844	1/ea, 10/bx, 20 bx/cs	A6212
7 in x 7 in, Sacral	886434	4845	1/ea, 10/bx, 20 bx/cs	A6213
9 in x 9 in, Sacral	1238306	4841	1/ea, 5/bx, 10 bx/cs	
9 in x 9 in, Heel	1264973	4839	1/ea, 5/bx, 10 bx/cs	



Hydrocellular Foam Dressing, Without Border, Silicone Adhesive, Sterile				
3 in x 3 in	1083085	4862	1/ea, 10/bx, 20 bx/cs	A6209
4 in x 4 in	1083086	4863	1/ea, 10/bx, 20 bx/cs	A6209
6 in x 6 in	1083087	4864	1/ea, 10/bx, 20 bx/cs	A6210
7 in x 7 in, Sacral	1083088	4865	1/ea, 10/bx, 20 bx/cs	A6211

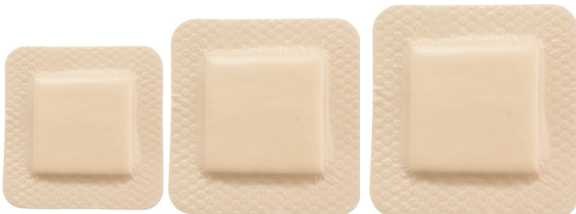

Lite Hydrocellular Foam Dressing, With Border, Silicone Adhesive, Sterile				
3 in x 3 in	1083089	4872	1/ea, 10/bx, 20 bx/cs	A6212
4 in x 4 in	1083090	4873	1/ea, 10/bx, 20 bx/cs	A6212
6 in x 6 in	1083091	4874	1/ea, 10/bx, 20 bx/cs	A6212



Lite Hydrocellular Foam Dressing, Without Border, Silicone Adhesive, Sterile				
3 in x 3 in	1083092	4892	1/ea, 10/bx, 20 bx/cs	A6209
4 in x 4 in	1083093	4893	1/ea, 10/bx, 20 bx/cs	A6209
6 in x 6 in	1083094	4894	1/ea, 10/bx, 20 bx/cs	A6210

Satisfaction Guaranteed
 For complete details, please visit
mms.mckesson.com/mckesson-brands.

Never Compromise on Care™

With McKesson Brands, you'll never have to compromise. We deliver quality products at competitive prices, backed by our 100% satisfaction guarantee. We have decades of experience in the market, and we offer you support and work to help ensure you have the products you need, when you need them—making it even easier for you to provide the best care for your patients.

The product information contained in this document, including the product images and additional product materials, was collected from various supplier sources. All product claims and specifications are those of the product suppliers and have not been independently verified by McKesson Medical-Surgical or its affiliates ("McKesson"). McKesson is not responsible for errors or omissions in the product information.

The properties of a product may change or be inaccurate following the posting or printing of the product information in the document, either in the print or online version. Caution should be exercised when using or purchasing any products from McKesson's online or print documents by closely examining the product packaging and the labeling prior to use.

Due to product changes, information listed in this document is subject to change without notice. This information is placed solely for your convenience in ordering and McKesson disclaims all responsibility for its completeness and accuracy, whether or not the inaccuracy or incompleteness is due to fault or error by McKesson.