Century Gold

Alternating and
Low Air Loss Pressure Relief System

User Manual



Manufactured by:
Caremed Supply, Inc.

7F., No.2, Ln. 235, Baoqiao Rd.,Xindian Dist.,
New Taipei City 231, Taiwan (R.O.C.)

Tel: 886-2-29179808
Fax: 886-2-29186505

Distributed by:
Prius Healthcare USA

160 Scarlet Blvd.
Oldsmar, FL 34677, USA

Tel: (813)854-5464
Fax: (813)854-5442

文件中心
2024.11.07
正式發行

IFUM29001 Ver A
2024/09/04


文件中心
2024.11.07
正式發行

Tabel of Contents

Safety Precautions 2
Other Safety Information 3
The Purpose of this Manual 4
Intended Use 4
Contraindications for use 4
Intended Users 4
Weight Capacity Guidelines 5
Key Features 5
Control Panel and Function 6
Specification 10
Symbols Used 12
Instructions for Use 13
Care and Cleaning 15
Waste Disposal 16
Storage and Care 17
Maintenance & Troubleshooting 18
Electromagnetic Compatibility 19
Expected Service Life 25
Warranty 25

Please finish reading this manual before using this equipment and keep this instruction handy.

文件中心
2024.11.07
正式發行

文件中心
2024.11.07
正式發行

Safety Precautions

When using electrical products, basic safety precautions should always be followed, including the following:

- This device can be used in home healthcare and professional healthcare environment.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State.

⚠ **CAUTION - READ ALL INSTRUCTIONS BEFORE USING THE APPLIANCE**

⚠ **DANGER** – To reduce the risk of electrocution:

1. Always unplug this device immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for the device that has fallen into water. Unplug immediately.
6. Do not use when the device is under repair or maintenance.

⚠ **WARNING** – To reduce the risk of burns, electrocution, fire or injury to persons:

1. Connect the Master Control unit to a proper power source. The device should never be left unattended when plugged in.
2. Close supervision is necessary when the device is used by, on, near children or physically challenged individuals.
3. Use the device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this device if it has a damaged cord or plug, not working properly, having been dropped or damaged, or dropped into water. Return the device to a service center for examination and repair.
5. Keep the cord away from heated surfaces and anywhere hard to unplug with the device like on the shelves, in a wardrobe or the socket on the ceiling.
6. Never block the air opening of this device or place it on a soft surface such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and other similar debris.
7. Never drop or insert objects into any openings.
8. Do not use this device outdoors.
9. Do not operate this device where aerosol (spray) products are being used or where oxygen / any flammable gases (such as Anesthetic Agents) are being administered.
10. DISCONNECT POWER SUPPLY BEFORE OPENING.
11. The device should only be repaired by authorized and qualified professional.
12. The product has no user/operator serviceable parts.

13. Do not modify the device without the authorization of the manufacturer. Failure to do so may exposure to hazards.
14. Maintain the device regularly. When the device needs repair or maintenance, please turn off the power first and unplug it from the mains, keep empty mattress, and do not use the device.
15. Keep the pump and mattress away from sources of liquid and open flames.
16. Keep the pump and mattress away from sharp objects.
17. If pain, irritation, numbness, swelling, or redness occurs discontinue use and contact a healthcare professional.
18. During DIP or power interruption the device would stop functioning and fail the intended use. Please ensure the device is used with a stable power /Uninterruptible Power Supply to prevent this scenario. Do not obstruct the mains plug or position the equipment where the connection to the mains line can be accidentally disconnected.
19. Power cable & pump shall be placed at the foot-side of the patient to prevent any risk of strangulation due to cable.
20. Medical electrical equipment needs special precautions regarding EMC and needs to be installed according to the EMC information provided.
21. Careful consideration of the EMC information is essential when stacking or collocating equipment and when routing cables and accessories.
22. To avoid the risk of cross-contamination, use liner or bed sheet on the mattress.
23. To avoid the risk of death or injury from falling, **side rails need to be in the raised** or guard position when a patient is on the bed.

⚠ **CAUTION** – When choosing a mattress, make sure the distance between top of the mattress (without compression) and top of the side rails and/or restraints is at least **8.66 in (220mm)** to help prevent inadvertent bed exit or falls.

Other Safety Information

Patient Entrance / Exit – It is essential for caregivers to assist patients in exiting the bed.

Brakes – Ensure that caster brakes are locked once the bed is positioned. Always check if the wheels are securely locked before the patient transfers to or from the bed.

Bed Height – To reduce the risk of falls or injury, maintain the bed at the lowest position when the patient is unattended.

Bed Elevation – To prevent patient migration, keep the head of the bed as low as possible.

文件中心
2024.11.07
正式發行

The Purpose of this Manual


This operation manual is mainly focused on the set up, cleaning and routine maintenance of the Century Gold System. We recommend you keeping this manual handy to answer most of the question related to the system.

Intended Use

This device can be used in home healthcare and professional healthcare environment.

Century Gold System is a unique and innovative specialized mattress replacement unit. The system utilizes a special low air loss technology with a high flow rate that provides pressure management for pressure ulcers prevention. The alternating function provides active prevention for pressure relief, especially for those in acute care and long-term care settings. The Century Gold System offers “deep-cell therapy”; whereby, the cells completely collapse providing zero pressure at the point of deflation. The soft-firm adjustment allows the patient to adjust the firmness or softness of the surface for optimal comfort through 8 digital scales.

Unique on demand air loss technology increase the efficiency of exhausting heat and humidity from hip area.

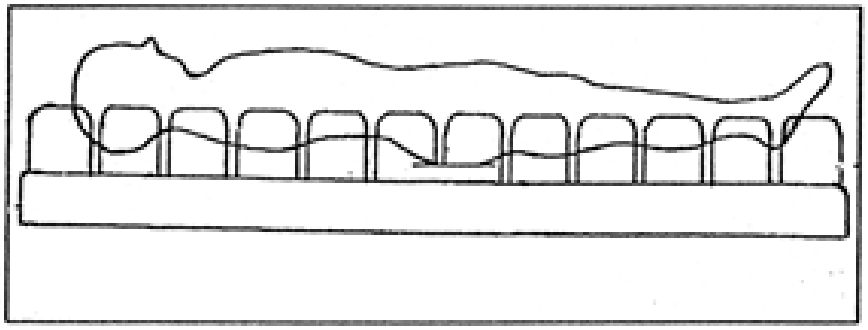

Contraindications for use

Alternating pressure therapy should not be used for patients with unstable fractures/spinal cord injuries, cervical or skeletal traction, gross oedema, burns or an intolerance to motion.

Intended Users

User(operator) of this device must have at least eight years of education and have the ability to read and understand English and Westernized Arabic numerals. No patient should be the operator of this device.

文件中心
2024.11.07
正式發行

Weight Capacity Guidelines

Mattress (applied part)			
Item No.	FM-PHR0057	FM-PHR0058	FM-PHR0060
Maximum weight	500 lbs (36”)	700 lbs (42”)	750 lbs (48”)

⚠ **CAUTION** – Do not exceed these weight limits at any time to ensure safety of the patient.


Key Features

- **Cycle Time** - alternating time can be adjusted in increments of 10, 15, 20 and 30 minutes.
- **Static Function** – during alternating function, the caregiver can select the “Static Function” to stop the alternation at any given time.
- **Side Rails** – the inflated side rails are designed to protect and prevent patients from falling out of bed. At the mid-section of the side rail, the top portion can be deflated for ease of patient exit.
- **Auto Firm** - the “Auto Firm” function provides uniform firmness for nursing procedures.
- **Audio & Visual Indicators** – adds protection to patient safety during power failures and low-pressure conditions. The audio indicator can be temporary turned off by pressing the “MUTE” button on the control panel. It will go off again the issue persist after 20 minutes.
- **On Demand Low Air Loss (LAL)** – to increase patient comfort by providing air circulation within the mattress, removing patient bodily moisture and lowering patient environment temperature.
- **Auto-Lock** – the control panel automatically locks if left untouched for 3 minutes to prevent accidental changes to modes and settings.
- **Fowler** – increases mattress pressure to prevent patient from bottoming out in upright positions.
- **Vapor Permeable Cover** – easy to clean, breathable cover, reducing friction and shear whilst optimizing skin microclimate.
- **CPR valve** – allows quick deflation for cardio-pulmonary resuscitation in emergency situations.
- **Cable Management** – cords are retained within the replacement mattress so the risk of tripping and cord pinching is eliminated.


文件中心
2024.11.07
正式發行

Control Panel and Function


Control Panel

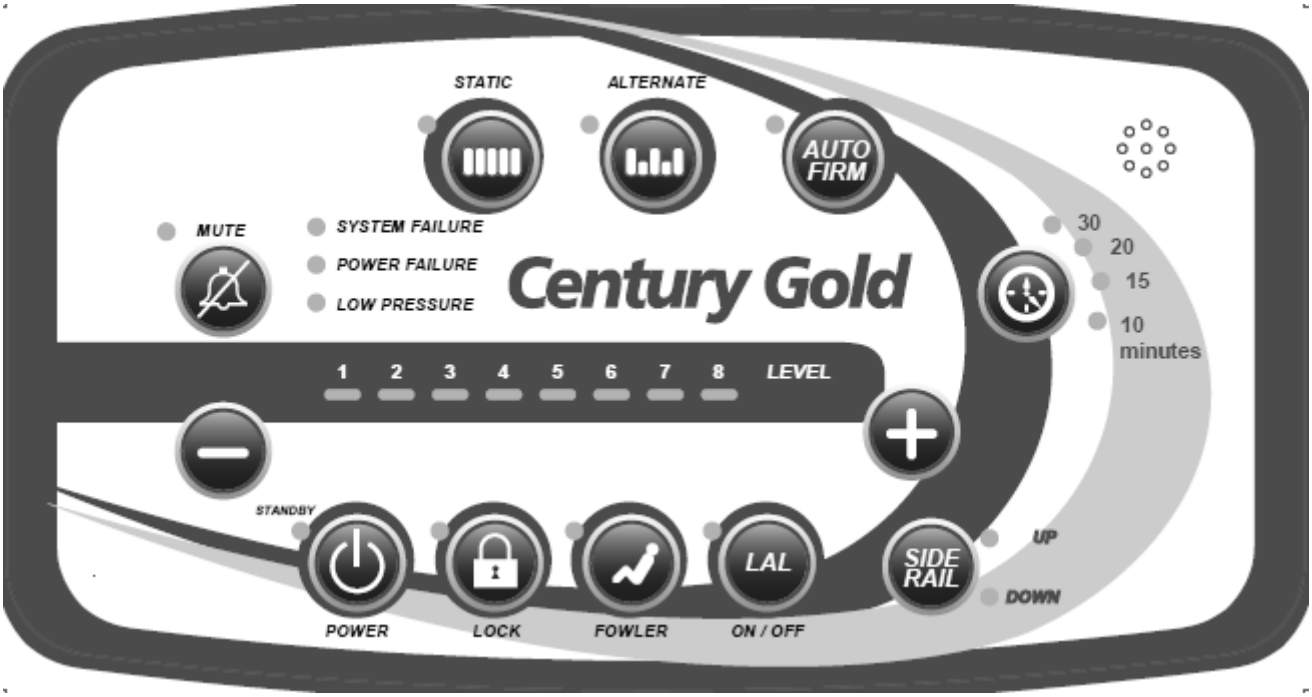


Functions:

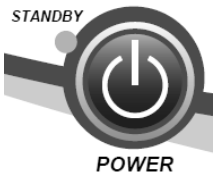	**Power On/off**: To turn on or off the control unit. After the device is plugged into a power source, the ‘STANDBY’ LED indicator will illuminate in amber. Press the ‘POWER’ button to activate the control unit; meanwhile, the ‘STANDBY’ LED indicator will be dimed.
	**Function Lockout:** Automatically locks out all functionality 3 minutes after a function change to ensure buttons cannot be inadvertently pressed. When lockout feature is activated, the lock LED indicator illuminate in green. Caregiver can unlock (or manually lock the device functions) by press and hold ‘LOCK’ button for 3 seconds.
	**Fowler:** When the patient is in upright position, the caregivers can press ‘FOWLER’ button to temporary increase the mattress pressure. The LED indicator will illuminate in green when activated. Press the ‘FOWLER’ button again to deactivate. The mattress will return to previous comfort setting. 

文件中心
2024.11.07
正式發行

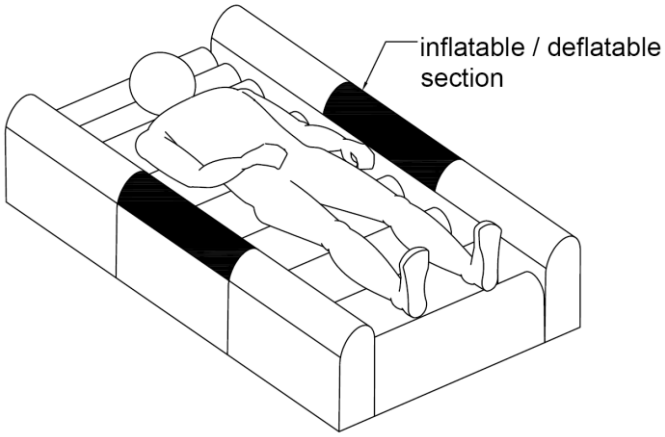	**On Demand Low Air Loss:** To increase patient comfort, the Low Air Loss function can be activated by pressing the 'LAL' button. The air will be circulated within the mattress to remove excessive moisture from the patient's body. When the function is in activate, the LED indicator will illuminate in green.
	**Side Rail:** The top portion of mid-section of the side rail can be deflated for ease of patient exit. Press the 'SIDE RAIL' button to inflate/deflate the side rail. When 'UP' LED indicator is in green, it means the rails are in raised position. Press the button again to lower the side rails. The 'DOWN' LED indicator will illuminate. **NOTE:** When switching between functions, the corresponding LED indicator will flash meaning the side rails is being inflated or deflated. It will become solid light when the inflation/deflation is complete. 
	**Comfort Level:** 8 comfort levels are available for a variety of patient weights. Use the (+) and (-) buttons to adjust pressure accordingly. When the device is initially activated, all comfort level LED indicators will flash in green signifying the mattress being inflated. This will continue for 20 mins to ensure the mattress is fully inflated, then it will return to default comfort level setting, and the LED indicator will illuminate in green.
	**Mute:** The control unit is equipped with both a visual and audio alert. When a failure is detected, the audio alert will sound. Press the 'MUTE' button to temporary silent the alert, meanwhile, the visual LED indicator will remain illuminated in amber. The alert will go off again if the problem persists in 20 minutes. Three types of failures can be detected: system failure, power failure, and low pressure.

文件中心
2024.11.07
正式發行

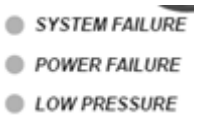

Failure modes: There are three types of failures that may be detected by the control unit, when a failure is detected, the failure mode LED indicator will illuminate in amber:

- System failure: abnormal pressure performance. Caregivers should turn off the controller and inspect all air cells, manifolds, and air hose for kinks.
- Power failure: check the power plug is securely plugged into a mains outlet.
- Low pressure: check the connections of the air hose and the control unit for proper connection.

NOTE: If the above measures are unsuccessful, contact your distributor for further instructions.




Static: Inflates all the air cells to the selected comfort level for ease of patient transfer or reposition. Press the 'STATIC' button to activate, and the LED indicator will illuminate in green.



Alternate: Press the 'ALTERNATE' button to enable 2-1 alternating mode, and the indicator will light up green. At least half of the cells will support patient weight at any time. Alternating time cycle can be adjusted, refer to 'Cycle Time' section (p. 9) for more details.

NOTE: It is important to perform hand check during alternating mode to avoid bottom out. Place the hand under the patient's buttock at the deflated cell section – the patient should have at least 5cm (3 fingers) of clearance between the buttocks and the bottom of the mattress. Finger check is recommended to be performed in every 2 hours.

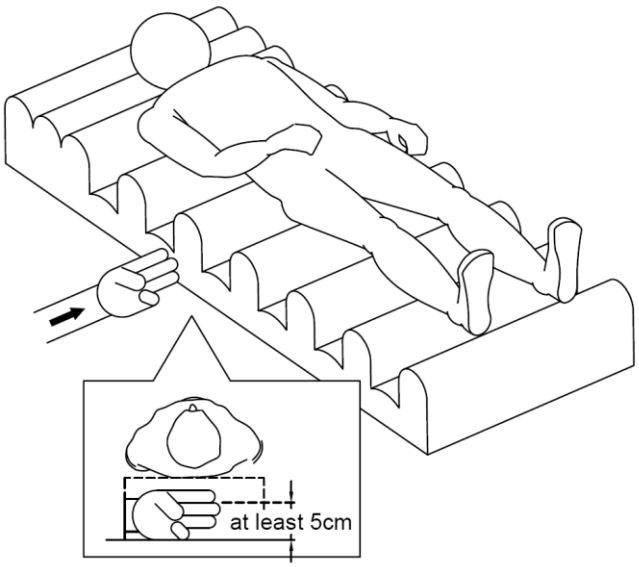

CAUTION:

Alternating pressure therapy is not recommended to patient who has serious pain or pain-sensitive symptom.

文件中心
2024.11.07
正式發行

	**Auto Firm**: Switches to static mode and Max-inflates all air cells that allows easy repositioning and transfers. This function can be initiated any time to fast inflate the mattress. Press the 'AUTO FIRM' button to activate the function, then the LED indicator will illuminate in green. Press the 'AUTO FIRM' button again to exit Auto Firm mode. For safety reason, the device will also automatically return to previous setting after 30 minutes. **NOTE:** When the control unit first initiated, it will automatedly enter into Auto Firm mode to fast inflate the mattress (the LED indicator will flash in green). Once it's fully inflated, it will return to static mode and set at pressure level 3.
	**Cycle time:** Cycle time is used only under alternation mode, and the default setting is 10 minutes. Press the 'CYCLE TIME' button to select the time setting (10, 15, 20, or 30 minutes) that's suitable for patient, and the corresponding LED indicator will illuminate in green.

CPR Valve: To quickly deflate the mattress in case of emergency. The CPR valve is located at patient's right hand side towards the head end. Turn the dial to 'OPEN' for quick deflation. Meanwhile, it is suggested to also disconnect the hose connector from the controller for air release at both head and foot end.

NOTE: The CPR valve should always stay closed under normal circumstance when the mattress is in use.

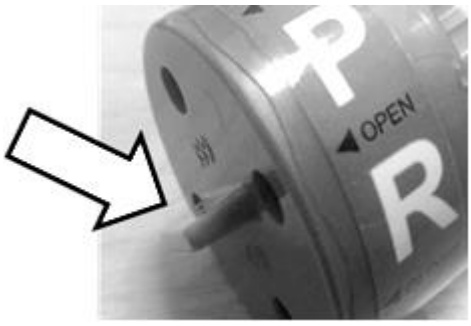


Fig. 1 Turn the dial to **OPEN** position to release air

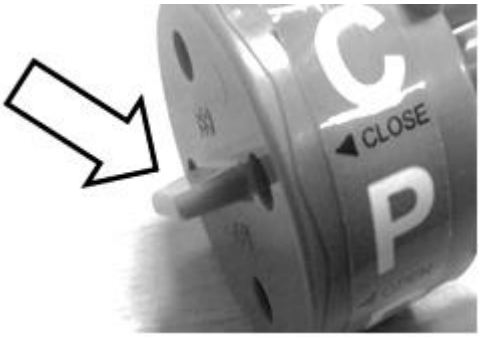


Fig. 2 Turn the dial to **CLOSE** position to inflate air mattress.

文件中心
2024.11.07
正式發行

Specification

Control Unit			
Model No:	M29-1		
Product Name	Century Gold Control Unit		
Item No.	FC-PHR0050		
Size in Inch (L x W x H):	15.7”x 8.27”x 10”		
Weight :	15 (lbs)		
Cycle Time:	10, 15, 20, 30 Minutes		
Max Operating Pressure	50mmHg		
Rated Voltage and Frequency	AC 100-240V, 50/60 Hz		
Max Current	0.8-0.4 A		
Fuse Rating	T2AH 250V		
Classification:	Class II, Type BF Not AP or APG type 		
Ingress of Water Protection:	IP21*		
Mode of Operation:	Continuous		
Power Cable	4.8m, non-shielding		
Environment (Temperature)	Operation: 15° C to 35° C (59° F to 95° F)		
	Storage and transportation: 5 ° C to 60 ° C (41° F to 140° F)		
Environment (Humidity)	Operation: 30% ~ 75% non-condensing		
	Storage and transportation: 30% ~ 90% non-condensing		
Operation Atmospheric Pressure Range	700 hPa to 1060 hPa		
Test Standard:	IEC60601-1, IEC60601-1-2 and IEC60601-1-11		
Mattress (applied part)			
Item No.	FM-PHR0057	FM-PHR0058	FM-PHR0060
Size in Inch (L x W x H)	84x36x10	84x42x10	84x48x10
Maximum weight	500 lbs (36”)	700 lbs (42”)	750 lbs (48”)
Cells Material	Nylon w/ PU		
Cover Type	Zipper cover		
Cover Material	Nylon woven fabric w/ PU coating finish		
Base Material	Woven Polyester fabric w/ PVC backing		

文件中心
2024.11.07
正式發行

*IP21

The IP Code (or International Protection Rating, sometimes also interpreted as Ingress Protection Rating) consists of the letters IP followed by two digits and an optional letter.

- First Digit: Solids

  The first digit indicates the level of protection that the enclosure provides against access to hazardous parts (e.g., electrical conductors, moving parts) and the ingress of solid foreign objects.



- Second Digit: Liquids

  Protection of the equipment inside the enclosure against harmful ingress of water.

IP Number	First Digit (Solids)	Second Digit (Liquids)
IP21	Protected from touched by fingers and objects greater than 12.5mm.	Against water: Vertical water drips.

文件中心
2024.11.07
正式發行

Symbols Used

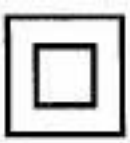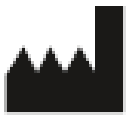	Type BF Protection Against Electronic Shock		Class II Equipment
	Refer To Instruction Manual / Booklet *NOTE: Background color: blue / Symbol: white*		Waste Disposal
	Caution, Consult accompanying documents		SGS product certification mark
	Date of manufacture		Manufacturer
	Batch code		Tumble Dry, Normal Low Heat
	Machine Wash. Water temperature should not exceed 40 ° C (104° F)		Do Not Bleach
	No Phenols		Do not dry clean
	Do Not Iron		

文件中心
2024.11.07
正式發行

Instructions for Use

WARNING – For important precautions, see p. 2~3.


WARNING – Check the power cord and plug for abrasions and excessive wear. Replace if necessary. Do NOT use the device if the power cord is not in a good condition.

CAUTION – Power on the device and activate ‘AUTO FIRM’. Check if the air hose and air cells have air leakage, or if air manifold is kinked or broken. Replace the parts if necessary.

CAUTION – Prior to usage, please allow the product to acclimate to ambient conditions for at least 30 minutes if the ambient temperature is 20°C. This helps ensure optimal performance and prevents potential damage caused by sudden temperature changes from the storage temperature range of 5°C to 60°C.

SETUP:

1. Replace the existing mattress with the Century Gold System. Orient the mattress so that the air hose is at the foot end of the bed.
2. Secure the straps beneath the mattress to the bed frame back board. Do NOT strap it to the side rail.

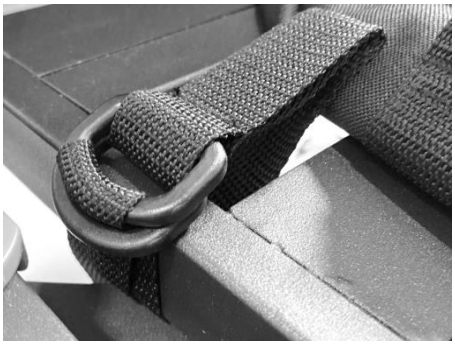3. Hang the control unit on the foot board of the bed frame, or place on the floor if there is no foot board attached to bed frame.
4. Attach the air hose to the socket located on the left side of the control unit. Make sure the connection is secure to avoid air leakage.
5. Check all air hoses to ensure no air hose is kinked.
6. Check if CPR valve is set to ‘CLOSE’ position. (Refer to CPR section, Fig. 2 in page 9)
7. Plug the power cable into the mains outlet, and the ‘STANDBY’ LED will illuminate amber.




8. Press the 'POWER' button on the control panel to inflate the mattress. While the mattress is being inflated, all the LED indicators for Auto Firm, Comfort Level, and Side Rail will flash in green.
9. It takes about 20 minutes for the mattress to be fully inflation. After the mattress is inflated, the system will enter into Static mode, where default comfort setting is set to level 3.
10. Caregivers now can transfer the patient onto the mattress.

文件中心
2024.11.07
正式發行

11. Please press and hold ‘LOCK’ button for 3 seconds to unlock the control functions.
12. Press the Side Rail button to ‘UP’ mode to prevent from patient falling. Then choose the suitable comfort level for the patient by pressing ‘SOFT’ (-) or ‘FIRM’ (+) button.
13. Press ‘ALTERNATE’ to activate 2-1 alternating mode. Perform finger check to ensure the distance between the patient’s buttock area and the bottom of the mattress is at least 5cm (about 3 fingers) – adjust the comfort level as needed. Finger check is recommended every 2 hours to prevent bottoming out.
14. Press the ‘CYCLE TIME’ button to select the preferred time settings (10, 15, 20 & 30 Minutes). Cycle time can only be use with Alternation mode, and the default setting is on 10 minutes.



CAUTION – Caregivers must check the controller setting and the patient condition every 2 hours. The bed frame side rails need to be in the upright position while the patient is on the mattress.

15. To stop the operation, press the Power’ button to turn off the control unit.
16. To pack the Century Gold System, follow the instruction guidelines as shown in Storage and Care section (p. 17).

文件中心
2024.11.07
正式發行

Care and Cleaning

WARNING – Turn off the control unit and disconnect the power plug from the mains outlet before cleaning. Failure to do so could result in personal injury or equipment damage.

WARNING – Do not expose the mattress or control unit to excessive moisture that would allow for liquid pooling. Personal injury or equipment damage could occur.

WARNING – Do not use harsh cleansers/detergents, such as sourcing pads and heavy-duty grease removers, or solvents, such as acetone. Equipment damage could occur.

WARNING – Do not use high temperature autoclave steam cleaning devices or phenolic based products for cleaning. Equipment damage could occur.

Cleaning and Infection Control

- To prevent from cross contamination, the mattress should be examined and disinfected between patient use.
- If the system is in constant use, it is recommended that the system is cleaned on a weekly basis.

Mattress

For basic care and cleaning, wipe down with damp soft cloth and mild detergent. The mattress top/bottom cover and mattress air cells can be wiped over with a solution of sodium hypochlorite 1000 ppm or any other non-phenolic germicidal solution. Dry thoroughly before reassembly.

If the top cover becomes severely soiled, remove the top cover from the bottom cover before washing. It can also be washed and thermally disinfected in a washing machine at a maximum of 104°F (40°C). Please refer to the top cover washing label for detailed cleaning instructions. Do NOT use the system without top cover.

Do NOT machine wash or tumble dry the air cells or bottom cover.

Control Unit

Disconnect control unit from mains power before cleaning.

Do not spray disinfectant directly onto the control unit, or immerse the unit in water or other fluid.

Gently wipe down the external case with a soft cloth using a solution of sodium hypochlorite 1000 ppm or any EPA-approved hospital-grade disinfectant. **Do NOT use phenol-based cleaning solutions for cleaning.**

文件中心
2024.11.07
正式發行

Air Filter

The air filter should be checked and changed regularly, at a minimum of every six months.

To replace the air filter, remove the filter cover on the back of the control unit, pinch the filter and pull it out from the back of the control unit. Place a new filter into the slot and secure the filter cover back onto the control unit.

Waste Disposal

This device has been supplied from an environmentally aware manufacturer that complies with the WEEE.

This device may contain substances that could be harmful to the environment if disposed of in places (landfills) that are not appropriate according the legislation. Please be environmentally responsible and recycle this product through your recycling facility at its end of life.



文件中心
2024.11.07
正式發行

Storage and Care

- Clean the Century Gold System prior to storage and prior to re-use.

To pack the Century Gold System, follow the instructions below:

1. Press the ‘Power’ button to turn off the control unit.
2. Disconnect the mattress hose connector from the control unit.
3. Unplug the power cord from the mains outlet, and then tighten it up with a cable strap.
4. Place the control unit inside a plastic bag, then back to the original carton.
5. To pack the mattress, turn the CPR valve to ‘OPEN’ position to expedite the air expulsion process (refer to CPR section, Fig. 1, p. 10).
6. Roll the mattress starting from the head towards the foot end, pausing to squeeze out any trapped air with each roll until it’s fully deflated.
7. Use the base mount straps to secure the rolled mattress. It can now be placed inside a plastic bag and back to the original carton.

It is recommended that the following guidelines are used whenever the system is being stored or transported to another location:

Temperature limitations	41°F~140°F (5°C ~ 60°C)
Relative humidity	30% ~ 90% non-condensing

文件中心
2024.11.07
正式發行

Maintenance & Troubleshooting

Daily maintenance is not required. Preventive maintenance is recommended to be carried out every 12 months. This equipment should only be serviced by authorized qualified technician. Parts and/or accessories supplied are specifically designed for use with this control unit. Use other products/parts in conjunction with the system is prohibited.

Troubleshooting:

Symptom	Inspection Procedure	Possible Solution
Air is pumping out from the control unit, but mattress is not inflating	Check if the connector of the air hose is connected properly	Re-connect the air hose
	Check if there is air leakage from the air hose	Replace with new hose.
	Check if there's any damage or air leakage on the air cells	Replace with new air cell(s)
	Check for kinks in air cell, air manifold, and air hose	Adjust the kinked part to enable smooth air flow
	Check if the CPR valve is left open	Close the CPR valve
The control unit is not functioning.	Check if the power cord is properly connected	Insert power plug to the mains outlet
	Check if the device is in standby mode	Turn the unit on by pressing POWER button
	Check the power cord for abrasions and excessive wear	Replace a new power cord.
	Check if the fuse is blown	Replace with a new fuse.
Some of the air cells are not properly inflated.	Check for kinks in air cells and/or air manifold	Adjust the kinked part to enable smooth air flow
	Check if there's any damage or air leakage on the air cells	Replace with new air cell(s)
No power failure visual and audio alert during power outage	When power resumes, check if the control unit is returned to normal operation. If so, leave the control unit plugged to a power source for 8 hours or more for battery charging. Then, unplug the control unit from the power source, and check if the 'POWER FAILURE' LED will flash in amber, accompanied by audio power failure alert.	If neither the visual nor audio alert activates, please contact your distributor for further instructions.

If the above measures are unsuccessful, contact your distributor for further instructions.

Warning: Do not attempt to repair the system yourself.

Warning: In case of damage, disconnect the plug from mains outlet. Do NOT use a damaged product.

文件中心
2024.11.07
正式發行

Electromagnetic Compatibility

WARNING: Medical electrical equipment needs special precautions regarding EMC and needs to be installed according to the EMC information provided. Careful consideration of this information is essential when stacking or collocating equipment and when routing cables and accessories.

WARNING: That portable RF communications equipment can affect medical electrical equipment. We recommend a safety distance no closer than 30 cm (12 inches) to any part of the M29-1 and at least 1 metre for sensitive equipment.

WARNING: Household electronic devices such as humidifiers, heaters or microwaves may cause interference with the device.

WARNING: X-rays, electric cauterizers, high-frequency, or high-energy devices may cause interference with the device.

CAUTION: If abnormal behavior is observed due to EM disturbances, please relocate the device accordingly.

CAUTION: Please do not use any other cables or accessories not approved by the manufacturer in this manual to avoid negative influence on electromagnetic compatibility.

Manufacturer's declaration-electromagnetic emissions		
The M29-1 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below. The customer or the user of the M29-1 should assure that it is used in such an environment.		
Emission test	**Compliance**	**Electromagnetic environment-guidance (for home and professional healthcare environment)**
RF emissions CISPR 11	Group 1	The M29-1 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The M29-1 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations /flicker emissions IEC 61000-3-3	Not applicable	

文件中心
2024.11.07
正式發行

Recommended separation distance between portable and mobile RF communications equipment and the M29-1

The M29-1 is intended for use in an electromagnetic environment (for home and professional healthcare) in which radiated RF disturbances are controlled. The customer or the user of the M29-1 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the M29-1 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	**Separation distance according to frequency of transmitter m**		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	**80 MHz to 800 MHz** $d = 1,2\sqrt{P}$	**800 MHz to 2,7 GHz** $d = 2,3\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance *d* in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

文件中心
2024.11.07
正式發行

Manufacturer's declaration-electromagnetic immunity

The M29-1 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.
The customer or the user of the M29-1 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance (for home and professional healthcare environment)
Electrostatic discharge(ESD) IEC 61000-4-2	Contact:±8 kV Air±2 kV,±4 kV,±8 kV,±15 kV	Contact:±8 kV Air±2 kV,±4 kV,±8 kV,±15 kV	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	± 2kV for power supply lines ± 1kV for input/output lines	± 2kV for power supply lines Not applicable	Mains power quality should be that of a typical home and professional healthcare environment.
Surge IEC 61000-4-5	± 0.5kV, ± 1kV line(s) to line(s) ± 0.5kV, ± 1kV, ± 2kV line(s) to earth	± 0.5kV, ± 1kV line(s) to line(s) Not applicable	Mains power quality should be that of a typical home and professional healthcare environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Voltage dips: 0 % *U*T; 0,5 cycle 0 % *U*T; 1 cycle 70 % *U*T; 25/30 cycles Voltage interruptions: 0 % *U*T; 250/300 cycle	Voltage dips: 0 % *U*T; 0,5 cycle 0 % *U*T; 1 cycle 70 % *U*T; 30 cycles Voltage interruptions: 0 % *U*T; 300 cycle	Mains power quality should be that of a typical home and professional healthcare environment. If the user of the M29-1 requires continued operation during power mains interruptions, it is recommended that the M29-1 be powered from an uninterruptible power supply.
Power frequency (50, 60 Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	30 A/m 60 Hz	The M29-1 power frequency magnetic fields should be at levels characteristic of a typical location in a typical home and professional healthcare environment.

NOTE UT is the a.c. mains voltage prior to application of the test level.

文件中心
2024.11.07
正式發行

Manufacturer's declaration-electromagnetic immunity

The M29-1 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.
The customer or the user of the M29-1 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance (for home and professional healthcare environment)
Conducted RF IEC 61000-4-6	3 Vrms: 0,15 MHz – 80 MHz 6 Vrms: in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz	3 Vrms: 0,15 MHz – 80 MHz 6 Vrms: in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz	**Portable and mobile RF communications equipment should be used no closer to any part of the** M29-1 including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	**Recommended separation distance:** $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P}$ 80MHz to 800 MHz $d = 2,3 \sqrt{P}$ 800MHz to 2,7 GHz Where *P* is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and *d* is the recommended separation distance in metres (m). Interference may occur in the vicinity of equipment marked with the following symbol: 

NOTE1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

文件中心
2024.11.07
正式發行

Manufacturer's declaration-electromagnetic immunity

Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment

The M29-1 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.
The customer or the user of the M29-1 should assure that it is used in such an environment.

Test frequency (MHz)	**Band** [a] (MHz)	**Service** [a]	**Modulation** [b]	**Maximum power** (W)	**Distance** (m)	**IMMUNITY TEST LEVEL** (V/m)	**Compliance LEVEL** (V/m) (for home and professional healthcare)
385	380 –390	TETRA 400	Pulse modulation b) 18 Hz	1,8	0,3	27	27
450	430 – 470	GMRS 460, FRS 460	FM c) ±5 kHz deviation 1 kHz sine	2	0,3	28	28
710	704 – 787	LTE Band 13, 17	Pulse modulation b) 217 Hz	0,2	0,3	9	9
745							
780							
810	800 – 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation b) 18 Hz	2	0,3	28	28
870							
930							
1720	1,700 – 1,990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation b) 217 Hz	2	0,3	28	28
1845							
1970							
2450	2,400 – 2,570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation b) 217 Hz	2	0,3	28	28
5240	5,100 – 5,800	WLAN 802.11 a/n	Pulse modulation b) 217 Hz	0,2	0,3	9	9
5500							
5785							

NOTE If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

a) For some services, only the uplink frequencies are included.

b) The carrier shall be modulated using a 50 % duty cycle square wave signal.

c) As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

文件中心
2024.11.07
正式發行

Manufacturer's declaration-electromagnetic immunity

Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields

The M29-1 is intended for use in the electromagnetic environment (for home and professional healthcare) specified below.

The customer or the user of the M29-1 should assure that it is used in such an environment.

Frequencies	Test Level [A/m]	Modulation	Dwell time [s]	Compliance LEVEL [A/m] (for home and professional healthcare)
30 kHz (a)	8	CW	3	8
134,2 kHz	65	Pulse modulation (b) 2,1 kHz	3	65 (c)
13,56 MHz	7,5	Pulse modulation (b) 50 kHz	3	7,5 (c)

Note:

(a) This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME AND PROFESSIONAL HEALTHCARE ENVIRONMENT.

(b) The carrier shall be modulated using a 50 % duty cycle square wave signal.

(c) r.m.s., before modulation is applied.

文件中心
2024.11.07
正式發行

Expected Service Life

- For maintaining basic safety and essential performance in regards to EMC, the M29-1 has an expected service life of five years. To maintain the condition of the device, service the system regularly according to the schedule recommended by Caremed Supply, Inc.
- Medical electrical equipment needs special precautions regarding EMC. Shall the device be used within one-mile distance from AM, FM, or TV broadcast antennas, it needs to be installed according to the EMC information provided.
- Do NOT use unapproved accessories or attempt to modify, disassemble or otherwise misuse the device or any of its components.

Warranty

- Caremed Supply, Inc. guarantees that this equipment is free from defects in materials and workmanship. Our obligation under this warranty is limited to the repair of equipment returned to the place of purchase within 12 months of delivery date.
- We agree to service / adjust any equipment returned, and to replace or repair any part that is proven to be a warranty defect, at no charge.
- This warranty excludes equipment damage through shipping, tampering, improper maintenance, carelessness, accident, negligence or misuse, or products that have been altered, repaired or dismantled other than with the manufacture's written authorization and by its approved procedures and by properly qualified technicians.
- In no event shall Caremed Supply, Inc. be liable for any direct, indirect or consequential damages or losses resulting from the use of equipment.

文件中心
2024.11.07
正式發行

Manufactured by:

Caremed Supply, Inc.

7F., No.2, Ln. 235, Baoqiao Rd., Xindian Dist.,

New Taipei City 231, Taiwan (R.O.C.)

Tel: 886-2-29179808

Fax: 886-2-29186505

Distributed by:

Prius Healthcare USA

160 Scarlet Blvd., Oldsmar, FL 34677, USA

TEL: (813)854-5464

FAX: (813)854-5442

文件中心
2024.11.07
正式發行