



**Whitehall Manufacturing®**  
Manufacturer of Healthcare and Rehabilitation Products since 1946

Division of  
Acorn Engineering Company®



A MEMBER OF



## Hi-Boy Whirlpools

Use



# H Series 60 Gallon Mobile Whirlpool

**WHH-60-M**



WHH-60-M

### Whitehall WHH-60-M Mobile Whirlpools

**Tank Capacity** is 60 gallons.

**Whirlpool Tank** is 36"L x 20"W x 28"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning.

**Heavy Duty Swivel Casters** support all mobile models. Rear pair of casters are lockable for stabilizing Whirlpool during use.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body during Whirlpool use. It is spring balanced and provided with a locking device.

**Turbine Assembly** functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

**Drain option -SDP** has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

**Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

### GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 60-gallon capacity, 36"L x 20"W x 28"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit [www.whitehallmfg.com](http://www.whitehallmfg.com)  
for most current specifications.



Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681

**WHH-60-M**

NEW: 09/10/2024



**Whitehall Manufacturing®**  
Manufacturer of Healthcare and Rehabilitation Products since 1946

## Hi-Boy Whirlpools

### MODEL NUMBER AND OPTION SELECTIONS

#### Series

☒ WHH Hi-Boy Series

#### Capacity

☒ -60 Gallons

#### Type

☒ -M Mobile Whirlpool

#### Turbine Electrical Requirements

☐ 115VAC / 50Hz ☐ 115VAC / 60Hz

☐ 230VAC / 50Hz ☐ 230VAC / 60Hz

Country of Use \_\_\_\_\_

#### For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet for details

#### Suffix Options (See Whirlpool Options Section for Details)



- ☐ -ADH Auxiliary Drain Hose Assembly
- ☐ -AS1 Adjustable Suspension Seat
- ☐ -EG Enviro-Glaze® Color Powder Coating - Specify Color: <https://www.whitehallmfg.com/enviro-glaze>
- ☐ -GFI Ground Fault Interrupter
- ☐ -LOGO Applied Logos: Team, Facility, Company Logo
- ☐ -SDP Separate Drain Pump
- ☐ -TTS1 Tank Top Seat

#### ELECTRIC REQUIREMENTS:

Amps: 6.9/6.2 Hertz: 50/60  
Volts: 115 GFI Receptacle

#### ELECTRIC REQUIREMENTS:

Amps: 3.5/3.1 Hertz: 50/60  
Volts: 230 GFI Receptacle

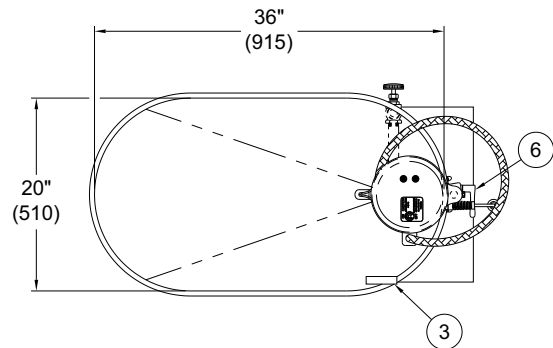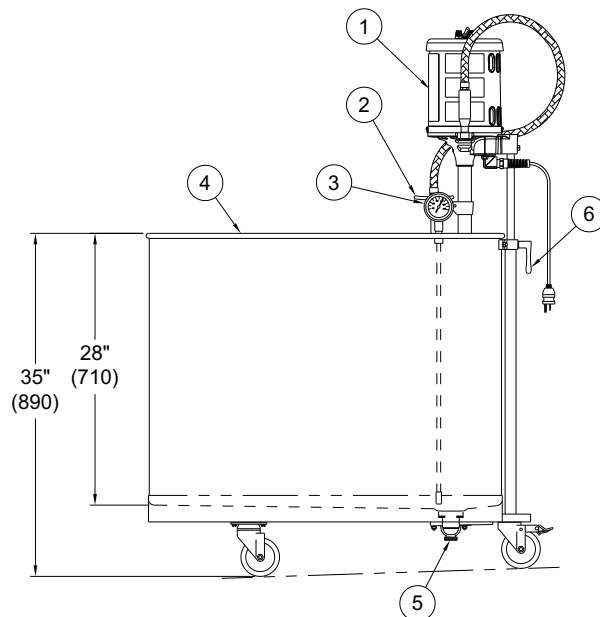

#### MOBILE WHH-60-M SERIES WHIRLPOOL

#### NOTES:

1. 1/2 hp Jet Pump Motor
2. Variable Pressure Control
3. Thermometer

4. Stainless Steel Tank
5. Auxiliary Drain Valve
6. Turbine Raising & Lowering Device

**WARNING:** Cancer and Reproductive Harm - [www.P65Warnings.ca.gov](http://www.P65Warnings.ca.gov)

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility for use of void or superseded data. © Whitehall Manufacturing, Member of Morris Group International. Please visit [www.whitehallmfg.com](http://www.whitehallmfg.com) for most current specifications.

#### SELECTION SUMMARY & APPROVAL FOR MANUFACTURING

Model Number & Options \_\_\_\_\_ Quantity \_\_\_\_\_  
Company \_\_\_\_\_ Date \_\_\_\_\_  
Contact \_\_\_\_\_ Title \_\_\_\_\_  
Approval for Manufacturing/Signature \_\_\_\_\_

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681

**WHH-60-M**

NEW: 09/10/2024