



Whitehall Manufacturing®
Manufacturer of Healthcare and Rehabilitation Products since 1946

Division of
Acorn Engineering Company®



A MEMBER OF



Lo-Boy Whirlpools

Use



L Series 105 Gallon Stationary Whirlpool with Legs

WHL-105-SL



WHL-105-SL

Whitehall WHL-105-SL Stationary Whirlpool with Legs

Whirlpool Tank is 66"L x 24"W x 18"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning. Whirlpool is fitted with adjustable stainless steel legs and a gravity drain valve.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body during Whirlpool use. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the

GUIDE SPECIFICATIONS

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 105-gallon capacity, 66"L x 24"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design. Whirlpool is fitted with adjustable stainless steel legs and a gravity drain valve. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com
for most current specifications.



Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681

WHL-105-SL

REV: 12/06/2024



Whitehall Manufacturing®
Manufacturer of Healthcare and Rehabilitation Products since 1946

Lo-Boy Whirlpools

MODEL NUMBER AND OPTION SELECTIONS

Series

☒ WHL Lo-Boy Series

Capacity

☒ -105 Gallons

Type

☒ -SL Stationary Whirlpool with Legs

Turbine Electrical Requirements

☐ 115VAC / 50Hz ☐ 115VAC / 60Hz

☐ 230VAC / 50Hz ☐ 230VAC / 60Hz

Country of Use _____

For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet for details.

Suffix Options (See Whirlpool Options Section for Details)


- ☐ -EG Enviro-Glaze® Color Powder Coating -
Specify Color: _____
<https://www.whitehallmfg.com/enviro-glaze>
- ☐ -GFI Ground Fault Interrupter
- ☐ -HR Head Rest, Adjustable
- ☐ -LOGO Applied Logos: Team, Facility, Company Logo
- ☐ -RS Reclining Seat
- ☐ -SDP Separate Drain Pump Assembly
- ☐ -TTS2 Tank Top Seat

ELECTRIC REQUIREMENTS:

Amps: 6.9/6.2 Hertz: 50/60
Volts: 115 GFI Receptacle

ELECTRIC REQUIREMENTS:

Amps: 3.5/3.1 Hertz: 50/60
Volts: 230 GFI Receptacle

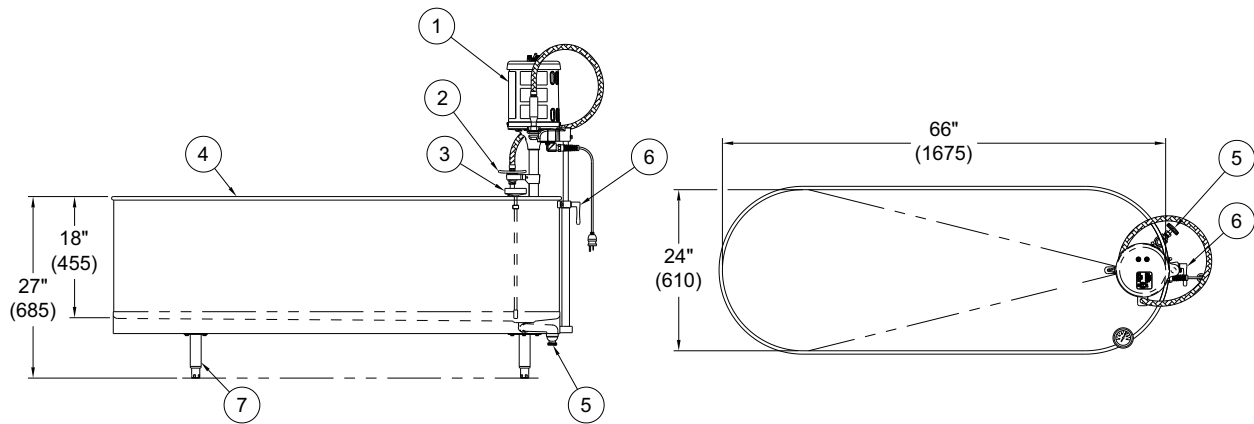

STATIONARY WHL-105-SL SERIES WHIRLPOOL

NOTES:

1. 1/2 hp Jet Pump Motor
2. Variable Pressure Control
3. Thermometer
4. Stainless Steel Tank

5. Auxiliary Drain Valve
6. Turbine Raising & Lowering Device
7. Adjustable Stainless Steel Legs

WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility for use of void or superseded data. © Whitehall Manufacturing, Member of Morris Group International. Please visit www.whitehallmfg.com for most current specifications.

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING

Model Number & Options _____ Quantity _____
Company _____ Date _____
Contact _____ Title _____
Approval for Manufacturing/Signature _____

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681

WHL-105-SL

REV: 12/06/2024